

Model B 28 | Standard-Incubators with mechanical adjustment

With the B series incubator, you get a reliable and powerful unit in a solid configuration at a very attractive price.

BENEFITS

- Safe and reproducible results even under high batch throughputs in long-term operation
- Inner glass door for a stable atmosphere



Model 28



Model 28

MAIN FEATURES

- Temperature range: 30 °C to 70 °C
- Hydro-mechanical thermostat
- Inner door made of tempered safety glass
- 2 chrome-plated racks
- stackable

ORDERING INFORMATION

Interior volume [cu.ft.]	Voltage	Option model	Version	Art.-No.
Model B 28				
1	230 V 1~ 50/60 Hz	Standard	B028-230V	9010-0002
		with TB class 1	B028-230V-T	9010-0004
	120 V 1~ 60 Hz	Standard	B028-120V	9010-0067

TECHNICAL DATA


Description	B028-230V ¹	B028-230V-T ¹	B028-120V ¹
Article Number	9010-0002	9010-0004	9010-0067
Performance Data Temperature			
Temperature range [°C]	30...70	30...70	30...70
Temperature variation at 37 °C [± K]	1.5	1.5	1.5
Temperature fluctuation at 37 °C [± K]	1	1	1
Heating-up time to 37 °C [min]	95	95	95
Electrical data			
Rated Voltage [V]	230	230	120
Power frequency [Hz]	50/60	50/60	60
Nominal power [kW]	0.25	0.25	0.25
Phase (Nominal voltage)	1~	1~	1~

¹ All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

Description	B028-230V ¹	B028-230V-T ¹	B028-120V ¹
Article Number	9010-0002	9010-0004	9010-0067
Outer dimensions			
Width net [in]	22.85	22.85	22.85
Height net [in]	15.96	15.96	15.96
Depth net [in]	16.74	16.74	16.74
Wall clearance back [in]	6.3	6.3	6.3
Wall clearance sidewise [in]	3.94	3.94	3.94
Doors			
Unit doors	1	1	1
Internal Dimensions			
Width [in]	15.76	15.76	15.76
Height [in]	11.03	11.03	11.03
Depth [in]	9.85	9.85	9.85
Measures			
Interior volume [cu.ft.]	1	1	1
Net weight of the unit (empty) [lbs]	51	51	51
permitted load [lbs]	55	55	55
Load per rack [lbs]	22	22	22
Fixtures			
Number of shelves (std./max.)	2/4	2/4	2/4

1 All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS INCL. FITTINGS AND CONNECTIONS [MM]

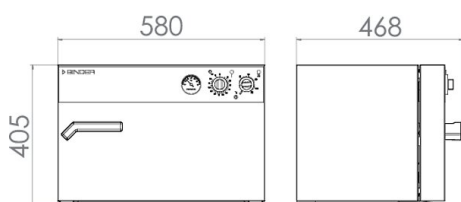

OPTIONS AND ACCESSORIES

Designation	Description	*	Art.-No.
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-0022
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	-	8012-0030
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	1002-0016
Rack	chrome plated	-	6004-0001
Rubber pads	set anti-slip feet	-	8012-0702
Shelf, perforated	Stainless steel	-	6004-0028

* Notes › See last page

SERVICES

Designation	Description	*	Art.-No.
Installation services			
Installation	and set up of unit at operating location, connect to existing connections	13, 18	DL10-0100
Instruction	unit function instructions for operation and programming of the controller	18	DL10-0500
Preventive maintenance			
Preventive maintenance	Executive of equipment inspection according to maintenance plan	14, 18	DL20-0200
Calibration services			
Calibration temperature incl. certificate	including certificate, one measuring point in center of chamber at specified temperature	14, 16, 17, 18	DL30-0101
Warranty service			
Extension of the warranty from 2 to 3 years	beginning with the date of delivery, wearing parts are not included	-	DL02-7041
Extension of the warranty from 2 to 5 years	beginning with the date of delivery, wearing parts are not included	-	DL02-7042

* Notes › See last page

NOTES

- 01 Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 02 UL mark is not granted when this option is used.
- 03 Heat resistant only to max. 200 °C.
- 04 Only available on units rated for 230 V.
- 06 Heating-up time may increase as a result of the lower heat conductivity.
- 07 The additional heat input may affect the temperature behavior.
- 09 Not in conjunction with the optional access port, door with window and interior lighting.
- 10 Not available on 23-liter units.
- 11 Not available on 23- or 53-liter units.
- 12 Only available on units rated for 230 V or 400 V.
- 13 Installation and connections take place at unit location; transport within the company only upon consultation.
- 14 We recommend a BINDER service contract to cover unit inspections, calibrations and validations.
- 15 OQ according to Yellow Paper = completed factory validation documentation of all OQ checklists.
- 16 Sensor calibration is performed in an accredited calibration laboratory.
- 17 Calibration is performed according to the BINDER factory standard.
- 18 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).
- 19 For additional accessories, refer to the Process documentation chapter.
- 20 When ordering IQ/OQ qualification folders and associated IQ/OQ execution on one order, we offer a 15 % discount for the item of the IQ/OQ folder.
When ordering IQ/OQ/PQ qualification folders and associated IQ/OQ/PQ execution on one order, we offer a 15 % discount for the item of the IQ/OQ/PQ folder.

BINDER GmbH
Tuttlingen, Germany
TEL +49 7462 2005 0
FAX +49 7462 2005 100
info@binder-world.com
www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.
Kowloon, Hong Kong, P.R. China
TEL +852 39070500
FAX +852 39070507
asia@binder-world.com
www.binder-world.com

Representative Office for CIS
Moscow, Russia
TEL +7 495 988 15 16
FAX +7 495 988 15 17
russia@binder-world.com
www.binder-world.ru

BINDER Environmental Testing
Equipment (Shanghai) Co., Ltd.
Shangha, P.R. China
TEL +86 21 685 808 25
FAX +86 21 685 808 29
china@binder-world.com
www.binder-world.com

BINDER Inc.
Bohemia, NY, USA
TEL +1 631 224 4340
FAX +1 631 224 4354
usa@binder-world.com
www.binder-world.us