

BINDER PURE AQUA SERVICE water treatment system



BINDER PURE AQUA SERVICE

KBF, KBF P, KBF LQC, KMF, KBWF, MKF, MKFT

The BINDER PURE AQUA SERVICE water treatment system is used for the treatment of tap water. It produces a water quality that is suitable for operating BINDER climate chambers. The BINDER PURE AQUA SERVICE water treatment system consists of a reusable measuring device for assessing water quality, a replaceable cartridge for water treatment, as well as the necessary connecting parts.

ADVANTAGES:

- Can be used irrespective of water quality
- It is not necessary to replace the resin; the cartridge is replaced when it is exhausted
- The shock-proof system is mobile and easy to transport
- Steady display of the current operating condition of the BINDER PURE AQUA SERVICE

TECHNICAL DATA

Inlet water temperature	5 – 40°C
Maximum flow rate (max.)	5 l/min
Supply pressure:	1 to 10 bar
Fresh water temperature:	+5°C...40°C
Permissible ingredients: max. permissible chlorine content:	0.1 mg/l


MAINS ADAPTER TECHNICAL DATA

Input	100 VAC to 240 VAC, 0.3 A, 50/60 Hz
Output	12 VDC, 500 mA
Protection class	IP 65 (dust-tight and protected against water jets).

SCOPE OF DELIVERY:

The BINDER PURE AQUA SERVICE water treatment system (item no. 8012-1809) consists of:

- PURE AQUA 300 single-use cartridge
- Mounting instructions
- Water hose, 3 m (metal hose with R3/4" connection) (OUT)
- Mains adapter, hose for drain cock, drain cock, seals conductivity measuring device

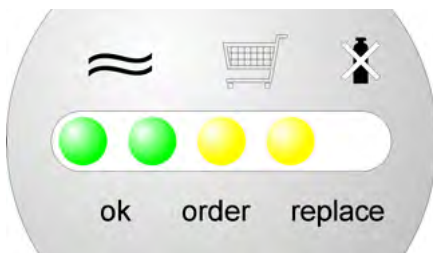

Conductivity measuring device with fill level indicator

BINDER PURE AQUA SERVICE COMPLETE SET	Height 385 mm, Ø 189 mm incl. PURE AQUA 300 single-use cartridge Item no. 8012-1809
PURE AQUA 300 SINGLE-USE REPLACEABLE CARTRIDGE	Height 385 mm, Ø 189 mm PURE AQUA 300 single-use cartridge Item no. 6011-0179
PURE AQUA 600 SINGLE-USE REPLACEABLE CARTRIDGE	Height 665 mm, Ø 189 mm PURE AQUA 600 single-use cartridge Item no. 6011-0180

► All extras online go2binder.com/en-Options