

# Model M 720 | Drying and heating chambers with forced convection and advanced program functions

With a maximum of 300 °C and various programming options the M series material test chambers are ideal for material and aging tests. With a high air flow, the powerful fans ensure rapid heating.

## BENEFITS

- Uniform drying conditions thanks to APT.line™ technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs
- Advanced programming options



Model 720



Model 720

## MAIN FEATURES

- Temperature range: +5 °C above ambient temperature to +300 °C
- APT.line™ preheating chamber technology
- Adjustable fan speed
- Program-controlled ventilation flap
- Color LCD controller with time-segment programming
- 2 chrome-plated racks
- Class 2 independent adjustable temperature safety device (DIN 12880) with visual alarm
- Computer interface: RS 422

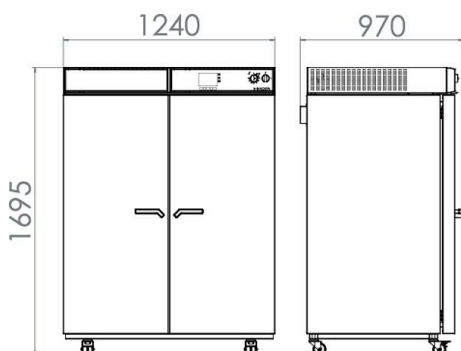
## ORDERING INFORMATION

Interior volume [L]	Voltage	Option model	Version	Art.-No.
720	400 V 3~ ph 50/60 Hz	Standard	M720-400V	9010-0205

**TECHNICAL DATA**

Description	M720-400V <sup>1</sup>
<b>Article Number</b>	<b>9010-0205</b>
<b>Performance Data Temperature</b>	
Temperature range +5 °C above ambient temperature to [°C]	300
Temperature fluctuation at 150 °C [± K]	0.3
Heating-up time to 150 °C [min]	21
Recovery time after 30 seconds door open at 150 °C [min]	3
<b>Air change data</b>	
Air change (approx.) at 150 °C [x/h]	36
<b>Electrical data</b>	
Rated Voltage [V]	400
Power frequency [Hz]	50/60
Nominal power [kW]	5
Phase (Nominal voltage)	3~
<b>Measures</b>	
Interior volume [L]	720
Net weight of the unit (empty) [kg]	203
Permitted load [kg]	120
Load per rack [kg]	45
Wall clearance back [mm]	160
Wall clearance sidewise [mm]	100
<b>Internal Dimensions</b>	
Width [mm]	1000
Height [mm]	1200
Depth [mm]	610
<b>Housing dimensions not incl. fittings and connections</b>	
Width net [mm]	1235
Height net [mm]	1690
Depth net [mm]	865
<b>Environment-specific data</b>	
Energy consumption at 150 °C [Wh/h]	1320
<b>Fixtures</b>	
Number of shelves (std./max.)	2/16

- 1 All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

**DIMENSIONS incl. fittings and connections [mm]**

**OPTIONS**

Designation	Description	*	Art.-No.
<b>Access port with silicone plug</b>	back		
	10 mm	01	8012-1287
	30 mm	01	8012-1371
	50 mm	01	8012-1387
	left		
	10 mm	01	8012-1280
	30 mm	01	8012-1363
	50 mm	01	8012-1379
	100 mm	01	8012-1511
	right		
	10 mm	01	8012-1273
	30 mm	01	8012-1358
	50 mm	01	8012-1375
	100 mm	01	8012-1508
	top		
	10 mm	01	8012-1266
30 mm	01	8012-1367	
50 mm	01	8012-1383	
<b>Analog output 4-20 mA</b>	for temperature values (output not adjustable)	-	8012-0480
<b>Calibration certificate, expanded</b>	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-1123
<b>Calibration certificate, temperature</b>	for temperature, measurement in center of chamber at specified temperature	-	8012-1142
	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-1601
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	-	8012-1580
	temperature measurement incl. certificate, 9 measuring points at specified temperature	-	8012-1559
<b>Door gasket</b>	made of FKM, silicone-free	-	8012-1315
<b>Door lock</b>	lockable door handle	-	8012-1256
<b>HEPA air filter</b>	On fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, > 99.995% at 0.3 µm)	-	8012-0222
<b>Inner chamber, reinforced</b>	including 2 reinforced racks max. total load 250 kg	-	8012-0343
<b>Lockable controller keyboard</b>	Lockable controller keyboard with key-operated switch	-	8012-0529
<b>Measurement of air exchange rate</b>	in accordance with ASTM D5374, definition and protocol according to ambient temperature	-	8012-1199
<b>Pt 100 temperature sensor</b>	additional flexible Pt 100, interior, for displaying the temperature on the unit display	-	8012-0334
<b>Viewing window and interior lighting</b>	viewing window in each door, 290 x 470 mm, and 30 W interior lighting 400 V option model	-	8012-0063

\* Notes › See last page

**ACCESSORIES**

Designation	Description	*	Art.-No.
<b>APT-COM™ 4 BASIC-Edition</b>	for simple logging and documentation requirements with up to 5 networked units. version 4, BASIC edition	-	9053-0039
<b>APT-COM™ 4 GLP-Edition</b>	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11. version 4, GLP edition	-	9053-0042
<b>APT-COM™ 4 PROFESSIONAL-Edition</b>	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units. version 4, PROFESSIONAL edition	-	9053-0040
<b>pH-neutral detergent</b>	concentrated, for gentle remove of residual contaminants; 1 kg	-	1002-0016

\* Notes › See last page

Designation	Description	*	Art.-No.
<b>Qualification documents</b>	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> , pressure, depending on unit		
	Digital in PDF format	–	7057-0001
	Hard copy inside folder	–	7007-0001
	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> – or pressure, depending on unit		
<b>Rack</b>	Digital in PDF format	–	7057-0005
	Hard copy inside folder	–	7007-0005
<b>Rack</b>	chrome plated	–	8012-2037
	stainless steel	–	8012-2082
<b>Rack accessories</b>	fasteners (1 set of 4) for additional security of racks	–	8012-0531
<b>Rack, reinforced</b>	stainless steel, with fasteners (1 set of 4)	–	8012-0374
<b>Shelf, perforated</b>	Stainless steel	–	8012-2081

\* Notes › See last page

## SERVICES

Designation	Description	*	Art.-No.
<b>Installation services</b>			
<b>Unit commissioning</b>	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	05	DL10-0100
<b>Unit instructions</b>	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05	DL10-0500
<b>Maintenance contracts</b>			
<b>Maintenance services</b>			
<b>BRONZE 3-year maintenance contract</b>	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts	05	DL20-0710
<b>GOLD 3-year maintenance contract</b>	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one test temperature specified by the user in the center of the usable space, including certificate	05	DL20-0910
<b>Maintenance</b>	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	05	DL20-0200
<b>SILVER 3-year maintenance contract</b>	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate	05	DL20-0810
<b>Calibration services</b>			
<b>Calibration certificate, temperature</b>	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate	03, 04, 05	DL30-0101
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate	03, 04, 05	DL30-0102
<b>Measurement of air exchange rate</b>	including certificate (in accordance with ASTM D5374)	03, 04, 05	DL33-0000
<b>Temperature measurement, 18 measuring points</b>	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-0118
<b>Temperature measurement, 27 measuring points</b>	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-0127
<b>Temperature measurement, 9 measuring points</b>	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-0109
<b>Validation services</b>			
<b>Execution of IQ/OQ</b>	Execution of IQ/OQ in accordance with qualification folder	05	DL41-0200
<b>Execution of IQ/OQ/PQ</b>	Execution of IQ/OQ/PQ in accordance with qualification folder	05	DL44-0500
<b>Warranty service</b>			
<b>1-year warranty extension</b>	The warranty is extended by 1 year from the delivery date, wear parts are excluded	–	DL50-0010

\* Notes › See last page

## NOTES

- 01 Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 02 UL mark is not granted when this option is used.
- 03 Sensor calibration is performed in an accredited calibration laboratory.
- 04 Calibration is performed according to the BINDER factory standard.
- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

### **BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
FAX +49 7462 2005 100  
info@binder-world.com  
www.binder-world.com

### **BINDER Asia Pacific (Hong Kong) Ltd.**

Kowloon, Hong Kong, P.R. China  
TEL +852 39070500  
FAX +852 39070507  
asia@binder-world.com  
www.binder-world.com

### **Representative Office for CIS**

Moscow, Russia  
TEL +7 495 988 15 16  
FAX +7 495 988 15 17  
russia@binder-world.com  
www.binder-world.com

### **BINDER Environmental Testing Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
FAX +86 21 685 808 29  
china@binder-world.com  
www.binder-world.com

### **BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
FAX +1 631 224 4354  
usa@binder-world.com  
www.binder-world.us