Model KBF-S 240 | Constant climate chambers with large temperature / humidity range

The Solid.Line KBF-S climate chamber is a custom-made solution for long-term testing and accelerated stability testing of pharmaceutical products in accordance with the ICH Q1A guideline. The corrosion-free and condensation-free interior also ensures exceptionally stable temperature and moisture levels, even when fully loaded.

BENEFITS

- Safe thanks to climatic homogeneity that far exceeds the accuracy required by ICH guidelines, even with a fully loaded unit.
- Reliable thanks to failsafe operation without compromise. The interior and Longlife evaporator plate are made entirely from stainless steel.
- · Smart, as a wide range of accessories makes it highly compatible for adaptation to specific customer requirements.
- Economical thanks to maximum effective volume which allows for 30% larger load compared to the competition.



Model 240



Model 240

MAIN FEATURES

- Temperature range: o °C to +70 °C
- Humidity range: 20% to 80% RH
- APT.line™ preheating chamber technology
- Humidity regulation with capacitive humidity sensor and steam humidification
- LCD to display temperature and humidity along with additional information and alarms
- Internal data logger, measured values can be read out in open format via USB
- Unit self-test for comprehensive status analysis
- Inner chamber made completely of stainless steel

- 2 stainless steel racks, incl. shelf supports
- Access port with silicone plug, 30 mm, left
- Robust twin casters with brakes, from 240 L
- Class 3.1 independent temperature safety device (DIN 12880) with visual and audible temperature alarm
- Computer interface: Ethernet
- Independent water supply via tank

ORDERING INFORMATION


| Interior volume [L] | Voltage | Option model | Version | ArtNo. |
|---------------------|-------------------------|--------------|----------------|-----------|
| 247 | 200230 V 1~ ph 50/60 Hz | Standard | KBFS240-230V | 9020-0366 |
| | 200240 V 1~ ph 50/60 Hz | Standard | KBFS240UL-240V | 9020-0367 |

TECHNICAL DATA


| Description | KBFS240-230V1 | KBFS240UL-240V1 |
|---|---------------|--|
| Article Number | 9020-0366 | 9020-0367 |
| Performance Data Temperature | | |
| Temperature range [°C] | 070 | 070 |
| Temperature variation at 40 °C [± K] | 0.3 | 0.3 |
| Temperature fluctuation [± K] | 0.1 | 0.1 |
| Max. heat compensation at 40 °C [W] | 300 | 300 |
| Performance Data Climate | | - |
| Temperature range [°C] | 1070 | 1070 |
| Temperature variation at 25 °C and 60 % RH [± K] | 0.3 | 0.3 |
| Temperature variation at 40 °C and 75 % RH [± K] | 0.3 | 0.3 |
| Temperature fluctuation at 25 °C and 60 % RH [± K] | 0.1 | 0.1 |
| Temperature fluctuation at 40 °C and 75 % RH [± K] | 0.1 | 0.1 |
| Humidity range | 2080 | 2080 |
| Humidity fluctuation at 25 °C and 60 % RH | 1,5 % RH | 1,5 % RH |
| Humidity fluctuation at 40 °C and 75 % RH | 1,5 % RH | 1,5 % RH |
| Recovery time after 30 seconds door open at 25 °C and 60 % RH [min] | 5 | 5 |
| Recovery time after 30 seconds door open at 40 °C and 75 % RH [min] | 11 | 11 |
| Electrical data | | |
| Rated Voltage [V] | 200230 | 200240 |
| Power frequency [Hz] | 50/60 | 50/60 |
| Nominal power [kW] | 1.4 | 1.4 |
| Unit fuse [A] | 16 | 16 |
| Phase (Nominal voltage) | 1~ | 1~ |
| Measures | • | <u>. </u> |
| Interior volume [L] | 247 | 247 |
| Net weight of the unit (empty) [kg] | 172 | 189 |
| Permitted load [kg] | 100 | 100 |
| Load per rack [kg] | 30 | 30 |
| Wall clearance back [mm] | 100 | 100 |
| Wall clearance sidewise [mm] | 100 | 100 |
| Internal Dimensions | | |
| Interior width [mm] | 650 | 650 |
| Interior height [mm] | 785 | 785 |
| Interior depth [mm] | 485 | 485 |
| Doors | 40) | 40) |
| Unit doors | 1 | 1 |
| Housing dimensions not incl. fittings and connections | • | - |
| Width net [mm] | 925 | 925 |
| Height net [mm] | 1460 | 1460 |
| Depth net [mm] | 800 | 800 |
| Environment-specific data | | |
| Energy consumption at 40 °C and 75 % RH [Wh/h] | 400 | 400 |
| Sound-pressure level [dB(A)] | 52 | 52 |
| Fixtures | <u>-ر</u> | <u>مر</u> |
| Number of shelves (std./max.) | 2/9 | 2/9 |
| number of stictics (stu./ max.) | 41 9 | 419 |

¹ All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS incl. fittings and connections [mm]

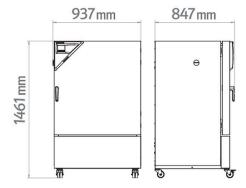

CHARTS

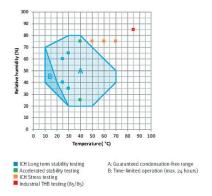

Climate chart

OPTIONS

| Designation | Description | * | ArtNo. |
|--------------------------------------|--|----|-----------|
| | left | | |
| | 30 mm | 01 | 8012-2072 |
| | 50 mm | 01 | 8012-2079 |
| | 100 mm | 01 | 8012-2077 |
| | right | | |
| access port with silicone | 30 mm | 01 | 8012-2071 |
| lug | 50 mm | 01 | 8012-2078 |
| | 100 mm | 01 | 8012-2076 |
| | top | | |
| | 30 mm | 01 | 8012-2073 |
| | 50 mm | 01 | 8012-2074 |
| | 100 mm | 01 | 8012-2075 |
| alibration certificate, expanded | for temperature and humidity; for extending the measurement in center of chamber to include another test value | = | 8012-2066 |
| | temperature measurement incl. certificate and 27 measuring points at specified temperature | - | 8012-2069 |
| Calibration certificate, temperature | temperature measurement incl. certificate, 15- 18 measuring points at specified temperature | - | 8012-2068 |
| | temperature measurement incl. certificate, 9 measuring points at specified temperature | - | 8012-2067 |
| alibration certificate, | Measurement in center of chamber at 25 °C / 60% RH or at specified test values | - | 8012-2065 |
| temperature and humidity | temperature (according to DIN1288o) and humidity measurement incl. certificate, 27 temperature measuring points and 1 humidity measuring point, at 25 $^{\circ}$ C / 60 $^{\circ}$ RH or at specified values | - | 8012-2070 |
| Door lock | lockable door handle | _ | 8012-2086 |

^{*} Notes > See last page

ACCESSORIES

| | · · · · · · · · · · · · · · · · · · · | | ArtNo. |
|--|--|---|-----------|
| APT-COM™ 4 BASIC- | for simple logging and documentation requirements with up to 5 networked units. | | |
| dition | version 4, BASIC edition | | 9053-0039 |
| APT-COM™ 4 GLP- Edition | for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11. | | |
| dition | version 4, GLP edition | - | 9053-004 |
| .РТ-СОМ™ 4 | convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units. | | |
| PROFESSIONAL-Edition | version 4, PROFESSIONAL edition | - | 9053-004 |
| BINDER PURE AQUA SERVICE | System for preparation or complete desalination of tap water, complete set containing PURE AQUA 300 single-use cartridge, measuring device, and all necessary connecting parts | | 8012-0759 |
| BINDER PURE AQUA SERVICE, accessories | Single-use, replacement cartridge for BINDER PURE AQUA System | | 6011-0165 |
| Convenient water supply | Flexible positioning of the water container with magnetic holder, hose adapter for water container with practical quick coupling, and water stop function | | |
| | External, with magnetic holder | - | 8012-1847 |
| | Basic set consisting of 2 pieces, attachment material, control unit for max. 4 light strips, 100-240 V, 50/60 Hz | | |
| | Basic set 300, length 30 cm | - | 8012-1107 |
| ED light bars | Basic set 500, length 50 cm | - | 8012-110 |
| ED light bars | Expansion set consisting of 2 pieces, attachment material: clips. For expanding the basic set of light bars | | |
| | Expansion set 300, length 30 cm | - | 8012-1716 |
| | Expansion set 500, length 50 cm | - | 8012-1717 |
| H-neutral detergent | concentrated, for gentle remove of residual contaminants; 1 kg | - | 8012-225 |
| | IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature and humidity values | | |
| | Digital in PDF format | - | 7057-000 |
| Qualification documents | Hard copy inside folder | - | 7007-000 |
| | IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature and humidity values | | |
| | Digital in PDF format | - | 7057-000 |
| | Hard copy inside folder | - | 7007-000 |
| tack | stainless steel | - | 8012-205 |
| tack accessories | fasteners (1 set of 4) for additional security of racks | - | 8012-228 |
| ack, reinforced | stainless steel, with fasteners (1 set of 4) | - | 8012-063 |
| helf, perforated | Stainless steel | - | 8012-216 |
| VLAN kit | The kit contains a client bridge that establishes a wireless connection between BINDER units and APT-COM4, LIMS, or customer-specific software via an Ethernet interface. It provides an alternative solution in situations where units need to be positioned in locations without an on-site Ethernet connection. In secured networks, installation and configuration must be performed by the customer's IT service. | - | 8012-226 |

Notes > See last page

SERVICES

| Designation | Description | * | ArtNo. |
|------------------------------------|--|----|-----------|
| Maintenance contracts | | | |
| Installation services | | | |
| Maintenance services | | | |
| Calibration services | | | |
| Validation services | | | |
| Warranty service | | | |
| 1-year warranty extension | The warranty is extended by 1 year from the delivery date, wear parts are excluded | - | DL50-0030 |
| BRONZE 3-year maintenance contract | Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts | 05 | DL20-0710 |
| Execution of IQ/OQ | Execution of IQ/OQ in accordance with qualification folder | 05 | DL42-0300 |
| Execution of IQ/OQ/PQ | Execution of IQ/OQ/PQ in accordance with qualification folder | 05 | DL44-0500 |
| | | | |

^{*} Notes > See last page

Data Sheet Model KBF-S 240

| Designation | Description | * | ArtNo. |
|--|--|---------------|-----------|
| GOLD 3-year maintenance contract | Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one temperature/humidity value, including certificate | 05 | DL20-0920 |
| Maintenance | One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate | 05 | DL20-0400 |
| SILVER 3-year maintenance contract | Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate | 05 | DL20-0820 |
| Temperature and humidity calibration | Expansion –Temperature and humidity calibration with 1 measuring point in center of chamber with 1 specified pair of values, including certificate | 03, 04, 05 | DL30-0302 |
| | Temperature and humidity calibration with 1 measuring point in center of chamber with 1 specified pair of values, including certificate | 03, 04, 05 | DL30-0301 |
| Temperature and humidity measurement according to DIN12880 | Temperature measurement in accordance with DIN 12880 with 27 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate | 03, 04, 05 | DL30-0427 |
| Temperature and humidity neasurement, 18-1 measuring points | Temperature measurement with 18 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate | 03, 04, 05 | DL30-0318 |
| Temperature and humidity neasurement, 27-1 measurement, 27-1 | Temperature measurement with 27 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate | 03, 04, 05 | DL30-0327 |
| Temperature and humidity neasurement, 9-1 measuring points | Temperature measurement with 9 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate | 03, 04, 05 | DL30-0309 |
| Jnit commissioning | Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) | 05 | DL10-0300 |
| Unit instructions | Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming | 05 | DL10-0700 |

^{*} Notes > See last page

NOTES

- Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge. UL mark is not granted when this option is used.
- 02
- Sensor calibration is performed in an accredited calibration laboratory. 03
- Calibration is performed according to the BINDER factory standard.
- Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany TEL +49 7462 2005 0 FAX +49 7462 2005 100 info@binder-world.com www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd. Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com

BINDER Inc. Bohemia, NY, USA TEL +1 631 224 4340 FAX +1 631 224 4354 usa@binder-world.com www.binder-world.us