


## Model BD 23 | Standard-Incubators with natural convection

The BINDER incubator of the BD series with natural convection is the specialist in long-term and stable continuous operation. This incubator is specially designed for gentle incubation of organisms, as well as conditioning of heat-sensitive media or microbiological heating.

### BENEFITS

- Safe and reproducible incubation even under high batch throughputs in long-term operation
- Disinfection routine at 100 °C
- Stable atmosphere by inner glass door

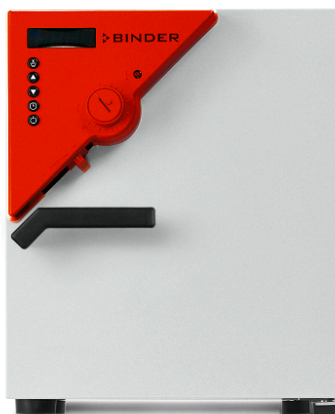

Model 23

### IMPORTANT FEATURES

- Temperature range: +5°C above ambient temperature to +100°C
- APT.line™ preheating chamber technology
- Natural convection
- Adjustable exhaust air flap
- Controller with timer functions
- Inner door made of tempered safety glass
- 2 chrome-plated racks
- Stackable
- Class 3.1 independent temperature safety device (DIN 12880) with visual temperature alarm
- Computer interface: RS 422

### ORDERING INFORMATION

Interior volume [L]	Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
<b>Model BD 23</b>					
20	230 V 1~ 50/60 Hz -12,5 A	CEE 7/7	Standard	BD023-230V	9010-0187
	120 V 1~ 60 Hz -12,5 A	NEMA 5-15	Standard	BD023UL-120V	9010-0189


## TECHNICAL DATA

Data		
Designation	BD023-230V	BD023UL-120V
Article Number	9010-0187	9010-0189
Option model	Standard	Standard
Performance Data Temperature		
Temperature range	+5 °C above ambient temperature to 100 °C	+5 °C above ambient temperature to 100 °C
Temperature uniformity at 37°C	0.5 ± K	0.5 ± K
Temperature fluctuation at 37°C	0.2 ± K	0.2 ± K
Heating up time to 37°C	59 min	59 min
Recovery time after door was opened for 30 s at 37°C	20 min	20 min
Electrical data		
Rated Voltage	230 V	120 V
Power frequency	50/60 Hz	60 Hz
Nominal power	0,2 kW	0,2 kW
Unit fuse	12,5 A	12,5 A
Phase (Nominal voltage)	1~	1~
Measures		
Interior volume	20 L	20 L
Net weight of the unit (empty)	27 kg	27 kg
Load per rack	12 kg	12 kg
Permitted load	25 kg	25 kg
Wall clearance back	100 mm	100 mm
Wall clearance sidewise	50 mm	50 mm
Housing dimensions not incl. fittings and connections		
Width net	435 mm	435 mm
Height net	495 mm	495 mm
Depth net	520 mm	520 mm
Internal Dimensions		
Interior width	222 mm	222 mm
Interior height	330 mm	330 mm

Interior depth	277 mm	277 mm
Unit doors	1	1
<b>Environment-specific data</b>		
Energy consumption at 37°C	11 Wh/h	11 Wh/h
<b>Fixtures</b>		
Number of shelves (std./max.)	2/4	2/4

All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

## DIMENSIONS Incl. fittings and connections [mm]

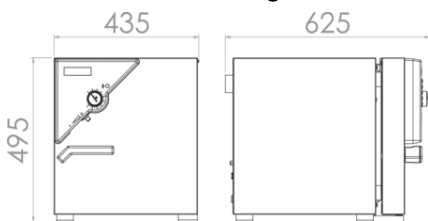

## OPTIONS

Designation	Description	BD 23	*	Art.-No.
	back			
	10 mm	•	01	8012-0390
	30 mm	•	01	8012-0397
	50 mm	•	01	8012-0399
	left			
	10 mm	•	01	8012-0388
	30 mm	•	01	8012-0047
Access port with silicone plug	50 mm	•	01	8012-0051
	right			
	10 mm	•	01	8012-0387
	30 mm	•	01	8012-0046
	50 mm	•	01	8012-0050
	top			
	10 mm	•	01	8012-0386
	30 mm	•	01	8012-0045

Designation	Description	BD 23	*	Art.-No.
	50 mm	•	01	8012-0049
	100 mm	•	01	8012-0053
Alarm function for overheating	Switchable acoustic alarm, with adjustable limit value on the independent temperature safety device			
	120 V option model	•	-	8012-0470
	230 V option model	•	-	8012-0469
Analog output 4-20 mA	for temperature values (output not adjustable)	•	-	8012-0432
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	-	8012-1110
	for temperature, measurement in center of chamber at specified temperature	•	-	8012-1129
Calibration certificate, temperature	temperature measurement incl. certificate, 9 measuring points at specified temperature	•	-	8012-0915
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	-	8012-0918
	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	-	8012-0921
Door lock	lockable door handle	•	-	8012-0318
Pt 100 temperature sensor	additional flexible Pt 100, interior, with external connection incl. LEMO connector (3-pin)	•	-	8012-0115

## ACCESSORIES

Designation	Description	BD 23	*	Art.-No.
APT-COM™ 4 BASIC-Edition	for simple logging and documentation requirements with up to 5 networked units. version 4, BASIC edition	•	-	9053-0039
APT-COM™ 4 GLP-Edition	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11. version 4, GLP edition	•	-	9053-0042
APT-COM™ 4 PROFESSIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units. version 4, PROFESSIONAL edition	•	-	9053-0040
Qualification documents	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> – or pressure, depending on unit			
	Digital in PDF format	•	-	7057-0005
	Hard copy inside folder	•	-	7007-0005

Designation	Description	BD 23	*	Art.-No.
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> , pressure, depending on unit			
	Digital in PDF format	•	-	7057-0001
	Hard copy inside folder	•	-	7007-0001
Rack	chrome plated	•	-	8012-2031
	stainless steel	•	-	8012-2046
Rubber pads	set anti-slip feet	•	-	8012-2030
Shelf, perforated	Stainless steel	•	-	8012-2164
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	•	-	9051-0018
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	-	8012-2250

## SERVICES

Designation	Description	*	Art.-No.
<b>Calibration services</b>			
Calibration certificate, temperature	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate	–	DL30-0101
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate	–	DL30-0102
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate	–	DL30-0109
Temperature measurement, 18 measuring points	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate	–	DL30-0118
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate	–	DL30-0127
<b>Installation services</b>			
Unit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	–	DL10-0100
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	–	DL10-0500
<b>Maintenance contracts</b>			
BRONZE 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts	–	DL20-0710
GOLD 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one test temperature specified by the user in the center of the usable space, including certificate	–	DL20-0910

Designation	Description	*	Art.-No.
SILVER 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate	–	DL20-0810
<b>Maintenance services</b>			
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	–	DL20-0200
<b>Validation services</b>			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	–	DL40-0100
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder	–	DL44-0500
<b>Warranty service</b>			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	–	DL50-0010

## NOTES

01      Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.

### **BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
FAX +49 7462 2005 100  
info@binder-world.com  
www.binder-world.com

### **BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
FAX +1 631 224 4354  
usa@binder-world.com  
www.binder-world.us

### **BINDER Environmental Testing Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
FAX +86 21 685 808 29  
china@binder-world.com  
www.binder-world.com

### **BINDER Asia Pacific (Hong Kong) Ltd.**

Kowloon, Hong Kong, P.R. China  
TEL +852 39070500  
FAX +852 39070507  
asia@binder-world.com  
www.binder-world.com