

Model KB PRO 1060 | Refrigerated incubators with an extended temperature range

A temperature control cabinet that can do it all, the KB PRO is our most powerful cooling incubator. With a temperature range of -20°C to +100°C and comprehensive programming options, it covers a wide range of applications that extend well beyond incubation. With optional light modules, the KB PRO can be converted into a plant growth chamber.

BENEFITS

- Reliable and powerful thanks to innovative temperature control technology
- Reliable: failsafe operation over the long term thanks to proven compressor technology and inner chamber made of stainless steel
- Smart: signal triangle with light-up status display, ergonomic control terminal with touchscreen program controller, rack with telescopic rails
- Efficient: low energy consumption thanks to modern inverter cooling technology



Model 1060



Model 1060

IMPORTANT FEATURES

- Temperature range: -20°C to +100°C
- Proven APT.line™ preheating chamber technology
- Energy-efficient cooling with modern inverter cooling unit
- Cooling with climate-neutral R-290 refrigerant
- Door heating to prevent condensation
- Adjustable fan speed
- Robust housing design with HIT sandwich technology
- Stainless-steel inner chamber
- Inner doors made of metal oxide-coated safety glass
- Two racks made of stainless steel (for models from 260 with telescopic rails)
- Access port Ø 30 mm with silicone plug (up to model 130)
- Access port 80 mm x 45 mm with silicone plug (from model 260)
- Tabletop units (up to model 130) are stackable with stacking adapter
- Floor-standing units (from model 260) featuring robust twin casters, with brakes at the front
- Intuitive touchscreen controller with time-segment and real-time programming
- Ergonomically adjustable control terminal (from model 470)
- Signal triangle with light-up status display (from model 470)
- Internal data logger, measured values can be read out in open format via USB
- Class 3.3 independent temperature safety device (DIN 12880) with visual and acoustic alarms
- Unit self-test for comprehensive status analysis
- Data interfaces: Ethernet, USB
- BINDER Multi Management Software APT-COM™ basic edition
- Optional light modules for ICH photostability testing (for models 260, 470, 720)
- Optional LED light modules for plant growth (for models from 260)

ORDERING INFORMATION

Interior volume		Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
[cu.ft.]	[L]					
Model KB PRO 1060						
38	1,070	208...240 V 1~ 50/60 Hz -16,0 A	CEE 7/7	Standard	KBPRO1060-230V	9020-0437


TECHNICAL DATA

Data	
Article Number	9020-0437
Designation	KBPRO1060-230V
Option model	Standard
Performance Data Temperature	
Temperature range	-20...100 °C
Temperature uniformity at 37°C	0.2 ± K
Temperature fluctuation at 37°C	0.1 ± K
Max. heat compensation at 37°C	1,000 W
Recovery time after door was opened for 30 s at 37°C	4 min
Electrical data	
Rated Voltage	208...240 V
Power frequency	50/60 Hz
Nominal power	1,8 kW
Unit fuse	16,0 A
Phase (Nominal voltage)	1~
Dimensions and weights	
Interior volume	38 cu.ft.
Net weight of the unit (empty)	778 lbs
Load per rack	88 lbs
Permitted load	441 lbs
Wall clearance back	3,94 in
Wall clearance sidewise	7,09 in
Housing dimensions not incl. fittings and connections	
Width net	50.59 in
Height net	76.81 in
Depth net	46.26 in
Internal Dimensions	
Interior width	39.8 in
Interior height	49.21 in

Interior depth	33.35 in
Inner doors	2
Unit doors	2
Environment-specific data	
Sound-pressure level at 25°C	50 dB(A)
Energy consumption at 25°C	325 Wh/h
Energy consumption at 37°C	340 Wh/h
Fixtures	
Number of shelves (std./max.)	2/16

All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS Incl. fittings and connections [mm]

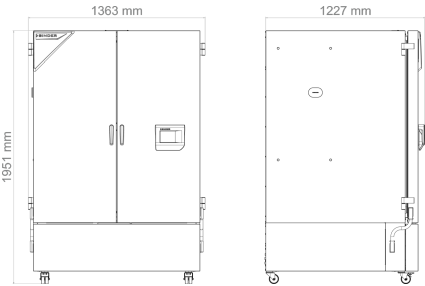

OPTIONS

Designation	Description	*	Art.-No.
Access port with silicone plug	left		
	30 mm	01	8012-2473
	50 mm	01	8012-2474
	100 mm	01	8012-2475
	right		
	30 mm	01	8012-2476
	50 mm	01	8012-2477
	100 mm	01	8012-2478
	top		
	30 mm	01	8012-2498

Designation	Description	*	Art.-No.
	50 mm	01	8012-2479
	100 mm	01	8012-2480
Alarm	accessible via 3-pin clamping sleeve (max. 24 V - 2.5 A)	-	8012-2502
Analog output 4-20 mA	for temperature values (output not adjustable)	-	8012-2573
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-2373
	for temperature, measurement in center of chamber at specified temperature	-	8012-2572
Calibration certificate, temperature	temperature measurement incl. certificate, 9 measuring points at specified temperature	-	8012-2377
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	-	8012-2381
	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-2425
Coated heat exchanger	For preventing damage to the aluminum heat exchanger caused by corrosive sample material	-	8012-2464
Door lock	lockable door handle	-	8012-2533
Door lock, electromechanical	can be opened with PIN code via controller	-	8012-2588
Inner chamber, reinforced	permitted load with U rails 300 kg	-	8012-2516
Interior lighting	With climate-resistant LEDs, switches when the unit door is opened	-	8012-2564
Shipping pallet	with unloading ramp for installing unit without forklift	-	8012-2546
Specimen temperature display	Flexible Pt 100 indoor temperature sensor, for object temperature display or control	-	8012-2553
Swivel casters, height-adjustable	For horizontal alignment of the unit on uneven floors	-	8012-2542

ACCESSORIES

Designation	Description	*	Art.-No.
APT-COM™ 4 GLP-Edition	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.		
	version 4, GLP edition	-	9053-0042
APT-COM™ 4 PROFESSIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.		
	version 4, PROFESSIONAL edition	-	9053-0040
Converter, voltage	to operate the device at a nominal voltage of 120 V 1 ph 60 Hz		
	240 V option model	-	8009-0821
CO ₂ gassing module	BINDER CO ₂ gassing module 0-1%, with NDIR CO ₂ sensor and analog output 4-20 mA	-	8012-2140

Designation	Description	*	Art.-No.
	BINDER CO ₂ gassing module 0-20%, with NDIR CO ₂ sensor and analog output 4-20 mA	-	8012-2138
LED light bars	Basic set consisting of 2 pieces, attachment material, control unit for max. 4 light strips, 100-240 V, 50/60 Hz		
	Basic set 300, length 30 cm	-	8012-1107
	Basic set 500, length 50 cm	-	8012-1108
	Basic set 900, length 90 cm	-	8012-1109
	Expansion set consisting of 2 pieces, attachment material: clips. For expanding the basic set of light bars		
	Expansion set 300, length 30 cm	-	8012-1716
	Expansion set 500, length 50 cm	-	8012-1717
	Expansion set 900, length 90 cm	-	8012-1718
Light module for plants	BINDER LED plant light module consisting of 16 LED strip lights, 3-channel controller and fresh air pump as well as independent power supply	-	8012-2439
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	8012-2250
Qualification documents IQ/OQ	IQ/OQ documents digital in PDF format – supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering:		
	temperature and light. Excluded: Document execution is separate	-	7057-0022
	IQ/OQ documents digital in PDF format – supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature. Excluded: Document execution is separate		
	temperature. Excluded: Document execution is separate	-	7057-0010
Qualification documents IQ/OQ/PQ	IQ/OQ/PQ documents digital in PDF format – supporting documents for validation customized to customer requirements, PQ section added to the qualification folder IQ/OQ; supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature values. Excluded: Document execution is separate		
	temperature values. Excluded: Document execution is separate	-	7057-0016
Reinforced rack	Shelf rack in set with pair of retaining rails (U rails), made entirely of stainless steel		
	Load-bearing capacity 80 kg	-	8012-2415
Shelf rack, perforated	Shelf rack in set with pair of retaining rails (U rails), made entirely of stainless steel		
	Load-bearing capacity 40 kg	-	8012-2409
Standard rack	Shelf rack in set with pair of retaining rails (ball bearing mounted telescopic rails), made entirely of stainless steel		
	Load-bearing capacity 40 kg	-	8012-2386
	Shelf rack in set with pair of retaining rails (U rails), made entirely of stainless steel		
	Load-bearing capacity 40 kg	-	8012-2392

SERVICES

Designation	Description	*	Art.-No.
Calibration services			
Calibration services	Calibration of CO ₂ concentration at 0,5 % by vol. with analyzed test gas, including certificate	05	DL30-0404
	Calibration of CO ₂ concentration at 5,0 % by vol. with analyzed test gas, including certificate	05	DL30-0405
Temperature calibration	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate Travel to be quoted separate.	05	DL30-0135
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate Travel to be quoted separate.	05	DL30-0102
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	05	DL30-0109
Temperature measurement, 18 measuring points	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 2 shelves are required.	05	DL30-0118
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	05	DL30-0127
Installation services			
Installation services	Installation of BINDER CO ₂ gassing module (0-1% or 0-20%), assembly of supply box, insertion of CO ₂ sensor and gas hose in interior, sealing of access port, connection of box to an existing CO ₂ gas supply, update of controller data set, function check	05	DL10-0130
	Installation of BINDER LED plant light module, assembly of supply box, assembly of LED strip lights as specified, laying of cables and fresh air hose, sealing of access port, function check	05	DL10-0150
Unit installation	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) Travel to be quoted as separate.	05	DL10-0170
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming. Travel to be quoted as separate Limited to 3 hours.	05	DL10-0550
Maintenance services			
BRONZE 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate. 10% Discount on Spare Parts and Labor for any additional Services needed Wear + tear parts to be purchased or billed to customer. Travel to be quoted as separate.	05	DL20-0720
GOLD 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate of Temperature. 25% Discount on Spare Parts and Labor for any additional Services needed. One additional courtesy visit per year can be used for service visit or operator. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	05	DL20-0945
Preventive maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Checking a test temperature specified by the user in center of usable	05	DL20-0604

Designation	Description	*	Art.-No.
	space without certificate. Wear + tear parts need to be purchased or billed to customer. Travel to be quoted as separate.		
SILVER 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate 15% Discount on Spare Parts and Labor for any additional Services needed One additional courtesy visit per year can be used for service visit or operator training. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	05	DL20-0845
Validation services			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder. Travel to be quoted separate.	05	DL41-0200
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder. Travel to be quoted separate.	05	DL44-0500
Warranty service			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	05	DL50-0010

NOTES

- 01 Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany
TEL +49 7462 2005 0
info@binder-world.com
www.binder-world.com

BINDER Inc.

Bohemia, NY, USA
TEL +1 631 224 4340
usa@binder-world.com
www.binder-world.us

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd.

Shanghai, P.R. China
TEL +86 21 685 808 25
china@binder-world.com
www.binder-world.com