

## Model KT 53 | Refrigerated incubators with Peltier technology

The KT series combines outstanding performance with impressive energy efficiency and environmental friendliness. The refrigerated incubators of the KT series are state-of-the-art in application and protect the samples.

### BENEFITS

- Safe and reproducible incubation conditions even at high environmental conditions
- Hot air disinfection at 100°C
- Energy efficient, environmentally friendly and secure



Model 53



Model 115

### IMPORTANT FEATURES

- Temperature range: 4 °C to 100 °C
- APT.line™ preheating chamber technology
- Adjustable fan speed
- Electric refrigeration with Peltier module
- Controller with time-segment and real-time programming
- Display via LCD monitor
- Input via pushbutton/rotary knob
- Inner door made of tempered safety glass
- 2 stainless steel racks
- Units up to 115 liters are stackable
- Class 3.1 independent temperature safety device (DIN 12880) with visual and audible temperature alarm
- Computer interface: Ethernet
- Data recording and USB interface

### ORDERING INFORMATION

| Interior volume |     | Power supply - unit fuse        | Plug*     | Version  | Model version | Art.-No.  |
|-----------------|-----|---------------------------------|-----------|----------|---------------|-----------|
| [cu.ft.]        | [L] |                                 |           |          |               |           |
| Model KT 53     |     |                                 |           |          |               |           |
| 2               | 53  | 200...240 V 1~ 50/60 Hz -10,0 A | CEE 7/7   | Standard | KT053-230V    | 9020-0311 |
|                 |     | 100...120 V 1~ 50/60 Hz -10,0 A | NEMA 5-20 | Standard | KT053UL-120V  | 9020-0312 |


## TECHNICAL DATA

| Data                                                  |             |              |
|-------------------------------------------------------|-------------|--------------|
| Article Number                                        | 9020-0311   | 9020-0312    |
| Designation                                           | KT053-230V  | KT053UL-120V |
| Option model                                          | Standard    | Standard     |
| Performance Data Temperature                          |             |              |
| Temperature range                                     | 4...100 °C  | 4...100 °C   |
| Temperature uniformity at 37°C                        | 0.3 ± K     | 0.3 ± K      |
| Temperature fluctuation at 37°C                       | 0.1 ± K     | 0.1 ± K      |
| Max. heat compensation at 40°C                        | 100 W       | 100 W        |
| Recovery time after door was opened for 30 s at 37°C  | 2 min       | 2 min        |
| Electrical data                                       |             |              |
| Rated Voltage                                         | 200...240 V | 100...120 V  |
| Power frequency                                       | 50/60 Hz    | 50/60 Hz     |
| Nominal power                                         | 0,4 kW      | 0,4 kW       |
| Unit fuse                                             | 10,0 A      | 10,0 A       |
| Phase (Nominal voltage)                               | 1~          | 1~           |
| Dimensions and weights                                |             |              |
| Interior volume                                       | 2 cu.ft.    | 2 cu.ft.     |
| Net weight of the unit (empty)                        | 139 lbs     | 139 lbs      |
| Load per rack                                         | 33 lbs      | 33 lbs       |
| Permitted load                                        | 88 lbs      | 88 lbs       |
| Wall clearance back                                   | 3,94 in     | 3,94 in      |
| Wall clearance sidewise                               | 9,45 in     | 9,45 in      |
| Housing dimensions not incl. fittings and connections |             |              |
| Width net                                             | 25.98 in    | 25.98 in     |
| Height net                                            | 25 in       | 25 in        |
| Depth net                                             | 24.8 in     | 24.8 in      |
| Internal Dimensions                                   |             |              |
| Interior width                                        | 15.75 in    | 15.75 in     |
| Interior height                                       | 15.75 in    | 15.75 in     |

|                               |          |          |
|-------------------------------|----------|----------|
| Interior depth                | 13.15 in | 13.15 in |
| Inner doors                   | 1        | 1        |
| Unit doors                    | 1        | 1        |
| Environment-specific data     |          |          |
| Sound-pressure level          | 48 dB(A) | 48 dB(A) |
| Energy consumption at 25°C    | 75 Wh/h  | 75 Wh/h  |
| Energy consumption at 37°C    | 75 Wh/h  | 75 Wh/h  |
| Fixtures                      |          |          |
| Number of shelves (std./max.) | 2/5      | 2/5      |

All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS Incl. fittings and connections [mm]

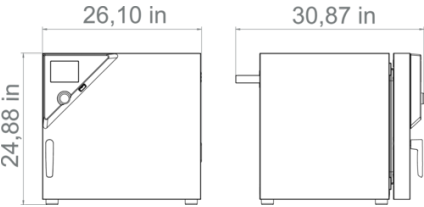

OPTIONS

| Designation                    | Description | KT 53 | *  | Art.-No.  |
|--------------------------------|-------------|-------|----|-----------|
| Access port with silicone plug | left        |       |    |           |
|                                | 10 mm       | •     | 01 | 8012-1290 |
|                                | 30 mm       | •     | 01 | 8012-1158 |
|                                | 50 mm       | •     | 01 | 8012-1165 |
|                                | right       |       |    |           |
|                                | 10 mm       | •     | 01 | 8012-1289 |
|                                | 30 mm       | •     | 01 | 8012-1155 |
|                                | 50 mm       | •     | 01 | 8012-1163 |
|                                | top         |       |    |           |
|                                | 10 mm       | •     | 01 | 8012-1288 |
|                                | 30 mm       | •     | 01 | 8012-1152 |
|                                | 50 mm       | •     | 01 | 8012-1161 |

| Designation                          | Description                                                                                                                                      | KT 53 | *  | Art.-No.  |
|--------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------|-------|----|-----------|
|                                      | 100 mm                                                                                                                                           | •     | 01 | 8012-1168 |
| Alarm output, zero-voltage           | accessible via 6-pin DIN socket (max. 24 V - 2.5 A)                                                                                              | •     | -  | 8012-0982 |
| Calibration certificate, expanded    | for temperature; for extending the measurement in center of chamber to include another test temperature                                          | •     | -  | 8012-1122 |
| Calibration certificate, temperature | for temperature, measurement in center of chamber at specified temperature                                                                       | •     | -  | 8012-1141 |
|                                      | temperature measurement incl. certificate, 9 measuring points at specified temperature                                                           | •     | -  | 8012-1558 |
|                                      | temperature measurement incl. certificate, 15- 18 measuring points at specified temperature                                                      | •     | -  | 8012-1579 |
|                                      | temperature measurement incl. certificate and 27 measuring points at specified temperature                                                       | •     | -  | 8012-1600 |
| Door lock                            | lockable door handle                                                                                                                             | •     | -  | 8012-1668 |
| Pt 100 temperature sensor            | additional flexible Pt 100, interior, for displaying the temperature on the unit display                                                         | •     | -  | 8012-0935 |
| Shelf, reinforced                    | positioned at bottom level, max. load 45 kg, with additional attachment for operation of shaking device, stirring device or roller bottle system | •     | -  | 8012-0818 |

## ACCESSORIES

| Designation                     | Description                                                                                                                                                                                                   | KT 53 | * | Art.-No.  |
|---------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|---|-----------|
| APT-COM™ 4 BASIC-Edition        | for simple logging and documentation requirements with up to 5 networked units.                                                                                                                               |       |   |           |
|                                 | version 4, BASIC edition                                                                                                                                                                                      | •     | - | 9053-0039 |
| APT-COM™ 4 GLP-Edition          | for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.                                                   |       |   |           |
|                                 | version 4, GLP edition                                                                                                                                                                                        | •     | - | 9053-0042 |
| APT-COM™ 4 PROFESSIONAL-Edition | convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.                                                                                                      |       |   |           |
|                                 | version 4, PROFESSIONAL edition                                                                                                                                                                               | •     | - | 9053-0040 |
| LED light bars                  | Basic set consisting of 2 pieces, attachment material, control unit for max. 4 light strips, 100-240 V, 50/60 Hz                                                                                              |       |   |           |
|                                 | Basic set 300, length 30 cm                                                                                                                                                                                   | •     | - | 8012-1107 |
|                                 | Expansion set consisting of 2 pieces, attachment material: clips. For expanding the basic set of light bars                                                                                                   |       |   |           |
|                                 | Expansion set 300, length 30 cm                                                                                                                                                                               | •     | - | 8012-1716 |
| pH-neutral detergent            | concentrated, for gentle remove of residual contaminants; 1 kg                                                                                                                                                | •     | - | 8012-2250 |
| Qualification documents IQ/OQ   | IQ/OQ documents digital in PDF format – supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature. Excluded: Document execution is separate |       |   |           |

| Designation                      | Description                                                                                                                                                                                                                                                                                                                                          | KT 53 | * | Art.-No.  |
|----------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|---|-----------|
|                                  | temperature. Excluded: Document execution is separate                                                                                                                                                                                                                                                                                                | •     | - | 7057-0010 |
| Qualification documents IQ/OQ/PQ | IQ/OQ/PQ documents digital in PDF format – supporting documents for validation customized to customer requirements, PQ section added to the qualification folder IQ/OQ; supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature values. Excluded: Document execution is separate |       |   |           |
|                                  | temperature values. Excluded: Document execution is separate                                                                                                                                                                                                                                                                                         | •     | - | 7057-0016 |
| Rack                             | stainless steel                                                                                                                                                                                                                                                                                                                                      | •     | - | 8012-2047 |
| Rack, reinforced                 | stainless steel, with fasteners (1 set of 4)                                                                                                                                                                                                                                                                                                         | •     | - | 8012-0829 |
| Rack accessories                 | fasteners (1 set of 4) for additional security of racks                                                                                                                                                                                                                                                                                              | •     | - | 8012-2280 |
| Rubber pads                      | set anti-slip feet                                                                                                                                                                                                                                                                                                                                   | •     | - | 8012-2030 |
| Shelf, perforated                | Stainless steel                                                                                                                                                                                                                                                                                                                                      | •     | - | 8012-2083 |
| Table on castors                 | stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm                                                                                                                                                                                                                                                                     | •     | - | 9051-0018 |

## SERVICES

| Designation                                  | Description                                                                                                                                                                                                                                                       | *  | Art.-No.  |
|----------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----|-----------|
| Calibration services                         |                                                                                                                                                                                                                                                                   |    |           |
| Temperature calibration                      | Calibration of one (1) test temperature specified by the user in center of chamber, including certificate Travel to be quoted separate.                                                                                                                           | 05 | DL30-0135 |
|                                              | Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate Travel to be quoted separate.                                                                                      | 05 | DL30-0102 |
| Temperature measurement, 9 measuring points  | Temperature measurement with 9 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.                                                                                   | 05 | DL30-0109 |
| Temperature measurement, 18 measuring points | Temperature measurement with 18 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 2 shelves are required.                                                                                  | 05 | DL30-0118 |
| Temperature measurement, 27 measuring points | Temperature measurement with 27 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.                                                                                  | 05 | DL30-0127 |
| Installation services                        |                                                                                                                                                                                                                                                                   |    |           |
| Unit installation                            | Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) Travel to be quoted as separate. | 05 | DL10-0170 |
| Unit instructions                            | Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming. Travel to be quoted as separate Limited to 3 hours.                                                                        | 05 | DL10-0550 |
| Maintenance services                         |                                                                                                                                                                                                                                                                   |    |           |

| Designation                        | Description                                                                                                                                                                                                                                                                                                                                                   | *  | Art.-No.  |
|------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----|-----------|
| BRONZE 3-year maintenance contract | Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate. 10% Discount on Spare Parts and Labor for any additional Services needed Wear + tear parts to be purchased or billed to customer. Travel to be quoted as separate.                                                                                                     | 05 | DL20-0720 |
| GOLD 3-year maintenance contract   | Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate of Temperature. 25% Discount on Spare Parts and Labor for any additional Services needed. One additional courtesy visit per year can be used for service visit or operator. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.      | 05 | DL20-0945 |
| Preventive maintenance             | One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Checking a test temperature specified by the user in center of usable space without certificate. Wear + tear parts need to be purchased or billed to customer. Travel to be quoted as separate. | 05 | DL20-0604 |
| SILVER 3-year maintenance contract | Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate 15% Discount on Spare Parts and Labor for any additional Services needed One additional courtesy visit per year can be used for service visit or operator training. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.              | 05 | DL20-0845 |
| Validation services                |                                                                                                                                                                                                                                                                                                                                                               |    |           |
| Execution of IQ/OQ                 | Execution of IQ/OQ in accordance with qualification folder. Travel to be quoted separate.                                                                                                                                                                                                                                                                     | 05 | DL41-0200 |
| Execution of IQ/OQ/PQ              | Execution of IQ/OQ/PQ in accordance with qualification folder. Travel to be quoted separate.                                                                                                                                                                                                                                                                  | 05 | DL44-0500 |
| Warranty service                   |                                                                                                                                                                                                                                                                                                                                                               |    |           |
| 1-year warranty extension          | The warranty is extended by 1 year from the delivery date, wear parts are excluded                                                                                                                                                                                                                                                                            | 05 | DL50-0010 |

**NOTES**

- 01      Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 05      Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

**BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
info@binder-world.com  
www.binder-world.com

**BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
usa@binder-world.com  
www.binder-world.us

**BINDER Environmental Testing  
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
china@binder-world.com  
www.binder-world.com