

Model KB 1060 | Refrigerated incubators with a large temperature range

The powerful virtuoso in cooling incubators for microorganisms: KB units are able to incubate with exceptional precision even when fully loaded, feature a hot-air disinfection routine, and can be programmed for weekend operation with ease. The new KB unit consumes as much as 40% less energy than the previous model.

BENEFITS

- Assured: highly precise incubation conditions thanks to power reserves
- Reliable: failsafe operation over the long term thanks to proven compressor technology and interior made of stainless steel
- Smart: hot-air disinfection and easy programming of a second set temperature
- Economical: up to 40% lower energy consumption than the previous model thanks to modern inverter cooling technology



Model 1060



Model 1060

IMPORTANT FEATURES

- Temperature range: 0 °C to +70 °C
- Proven APT.line™ preheating chamber technology
- Energy-efficient cooling with modern inverter cooling unit
- Cooling with climate-neutral R-600a refrigerant
- Door heating to prevent condensation
- Sturdy housing design with HIT sandwich technology
- Interior made completely of stainless steel
- Inner doors made of metal oxide-coated safety glass
- Two racks made of stainless steel
- 80 mm x 45 mm access port with silicone plug
- Tabletop units are stackable with stacking adapter
- Floor-standing units featuring sturdy twin casters with brakes at the front
- Simple programming of a second set temperature value
- Internal data logger, measured values can be read out in open format via USB
- Class 3.3 independent temperature safety device (DIN 12880) with visual and acoustic alarms
- Unit self-test for comprehensive status analysis
- Data interfaces: Ethernet, USB

ORDERING INFORMATION

Interior volume		Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
[cu.ft.]	[L]					
Model KB 1060						
38	1,070	220...240 V 1~ 50/60 Hz -16,0 A	CEE 7/7	Standard	KB1060-230V	9020-0476
		120 V 1~ 50/60 Hz -16,0 A	NEMA 5-20	Standard	KB1060UL-120V	9020-0504


TECHNICAL DATA

Data		
Article Number	9020-0476	9020-0504
Designation	KB1060-230V	KB1060UL-120V
Option model	Standard	Standard
Performance Data Temperature		
Temperature range	0...70 °C	0...70 °C
Temperature uniformity at 37°C	0.2 ± K	0.2 ± K
Temperature fluctuation at 37°C	0.1 ± K	0.1 ± K
Max. heat compensation at 37°C	640 W	640 W
Recovery time after door was opened for 30 s at 37°C	3 min	3 min
Electrical data		
Rated Voltage	220...240 V	120 V
Power frequency	50/60 Hz	50/60 Hz
Nominal power	1,5 kW	1,5 kW
Unit fuse	16,0 A	16,0 A
Phase (Nominal voltage)	1~	1~
Dimensions and weights		
Interior volume	38 cu.ft.	38 cu.ft.
Net weight of the unit (empty)	694 lbs	694 lbs
Load per rack	88 lbs	88 lbs
Permitted load	441 lbs	441 lbs
Wall clearance back	3,94 in	3,94 in
Wall clearance sidewise	7,09 in	7,09 in
Housing dimensions not incl. fittings and connections		
Width net	50.59 in	50.59 in
Height net	76.81 in	76.81 in
Depth net	46.26 in	46.26 in
Internal Dimensions		
Interior width	39.8 in	39.8 in
Interior height	49.21 in	49.21 in

Interior depth	33.35 in	33.35 in
Inner doors	2	2
Unit doors	2	2
Environment-specific data		
Sound-pressure level at 25°C	52 dB(A)	52 dB(A)
Energy consumption at 25°C	240 Wh/h	270 Wh/h
Energy consumption at 37°C	260 Wh/h	300 Wh/h
Fixtures		
Number of shelves (std./max.)	2/16	2/16

All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS Incl. fittings and connections [mm]

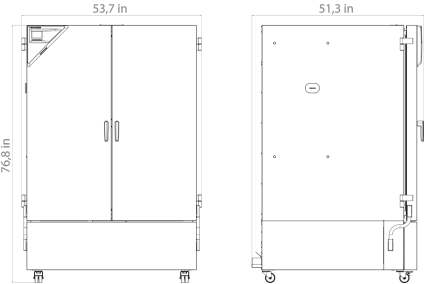

OPTIONS

Designation	Description	KB 1060	*	Art.-No.
Access port with silicone plug	left			
	30 mm	•	01	8012-2465
	50 mm	•	01	8012-2466
	100 mm	•	01	8012-2467
	right			
	30 mm	•	01	8012-2468
	50 mm	•	01	8012-2469
	100 mm	•	01	8012-2470
	top			
	30 mm	•	01	8012-2497

Designation	Description	KB 1060	*	Art.-No.
	50 mm	•	01	8012-2471
	100 mm	•	01	8012-2472
Alarm	accessible via 3-pin clamping sleeve (max. 24 V - 2.5 A)	•	-	8012-2370
Analog output 4-20 mA	for temperature values (output not adjustable)	•	-	8012-2371
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	-	8012-2372
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	•	-	8012-2571
	temperature measurement incl. certificate, 9 measuring points at specified temperature	•	-	8012-2376
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	-	8012-2380
	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	-	8012-2424
Door lock	lockable door handle	•	-	8012-2530
Door lock, electromechanical	can be opened with PIN code via controller	•	-	8012-2587
Inner chamber, reinforced	permitted load with U rails 300 kg	•	-	8012-2510
Interior lighting	With climate-resistant LEDs, switches when the unit door is opened	•	-	8012-2563
Shipping pallet	with unloading ramp for installing unit without forklift	•	-	8012-2545
Specimen temperature display	Flexible Pt 100 indoor temperature sensor, for object temperature display or control	•	-	8012-2549
Swivel casters, height-adjustable	For horizontal alignment of the unit on uneven floors	•	-	8012-2541

ACCESSORIES

Designation	Description	KB 1060	*	Art.-No.
APT-COM™ 4 BASIC-Edition	for simple logging and documentation requirements with up to 5 networked units.			
	version 4, BASIC edition	•	-	9053-0039
APT-COM™ 4 GLP-Edition	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.			
	version 4, GLP edition	•	-	9053-0042
APT-COM™ 4 PROFESSIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.			
	version 4, PROFESSIONAL edition	•	-	9053-0040
LED light bars	Basic set consisting of 2 pieces, attachment material, control unit for max. 4 light strips, 100-240 V, 50/60 Hz			

Designation	Description	KB 1060	*	Art.-No.
	Basic set 500, length 50 cm	•	-	8012-1108
	Basic set 900, length 90 cm	•	-	8012-1109
	Expansion set consisting of 2 pieces, attachment material: clips. For expanding the basic set of light bars			
	Expansion set 500, length 50 cm	•	-	8012-1717
	Expansion set 900, length 90 cm	•	-	8012-1718
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	-	8012-2250
Qualification documents IQ/OQ	IQ/OQ documents digital in PDF format – supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature. Excluded: Document execution is separate			
	temperature. Excluded: Document execution is separate	•	-	7057-0010
Qualification documents IQ/OQ/PQ	IQ/OQ/PQ documents digital in PDF format – supporting documents for validation customized to customer requirements, PQ section added to the qualification folder IQ/OQ; supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature values. Excluded: Document execution is separate			
	temperature values. Excluded: Document execution is separate	•	-	7057-0016
Reinforced rack	Shelf rack in set with pair of retaining rails (U rails), made entirely of stainless steel			
	Load-bearing capacity 80 kg	•	-	8012-2415
Shelf rack, perforated	Shelf rack in set with pair of retaining rails (U rails), made entirely of stainless steel			
	Load-bearing capacity 40 kg	•	-	8012-2409
Standard rack	Shelf rack in set with pair of retaining rails (ball bearing mounted telescopic rails), made entirely of stainless steel			
	Load-bearing capacity 40 kg	•	-	8012-2386
	Shelf rack in set with pair of retaining rails (U rails), made entirely of stainless steel			
	Load-bearing capacity 40 kg	•	-	8012-2392

SERVICES

Designation	Description	*	Art.-No.
Calibration services			
Temperature calibration	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate Travel to be quoted separate.	05	DL30-0135
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate Travel to be quoted separate.	05	DL30-0102
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	05	DL30-0109

Designation	Description	*	Art.-No.
Temperature measurement, 18 measuring points	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 2 shelves are required.	05	DL30-0118
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	05	DL30-0127
Installation services			
Unit installation	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) Travel to be quoted as separate.	05	DL10-0170
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming. Travel to be quoted as separate Limited to 3 hours.	05	DL10-0550
Maintenance services			
BRONZE 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate. 10% Discount on Spare Parts and Labor for any additional Services needed Wear + tear parts to be purchased or billed to customer. Travel to be quoted as separate.	05	DL20-0720
GOLD 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate of Temperature. 25% Discount on Spare Parts and Labor for any additional Services needed. One additional courtesy visit per year can be used for service visit or operator. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	05	DL20-0945
Preventive maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Checking a test temperature specified by the user in center of usable space without certificate. Wear + tear parts need to be purchased or billed to customer. Travel to be quoted as separate.	05	DL20-0604
SILVER 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate 15% Discount on Spare Parts and Labor for any additional Services needed One additional courtesy visit per year can be used for service visit or operator training. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	05	DL20-0845
Validation services			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder. Travel to be quoted separate.	05	DL41-0200
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder. Travel to be quoted separate.	05	DL44-0500
Warranty service			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	05	DL50-0010

NOTES

- 01 Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany
TEL +49 7462 2005 0
info@binder-world.com
www.binder-world.com

BINDER Inc.

Bohemia, NY, USA
TEL +1 631 224 4340
usa@binder-world.com
www.binder-world.us

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd.

Shanghai, P.R. China
TEL +86 21 685 808 25
china@binder-world.com
www.binder-world.com