

Model UF V 700 | Ultra-low temperature freezers with climate-neutral refrigerants


The BINDER ultra low temperature freezer ensures the safe storage of samples at -90 °C. It combines outstanding environmental friendliness with low energy consumption, convenient operation and an individual safety concept.

BENEFITS

- Lowest energy consumption in its class
- Environmentally-friendly refrigerants R-290 and R-170
- Thermal insulation with unique long-life VIPs (vacuum insulation panels)
- Rust-proof interior made completely of stainless steel



Model 700

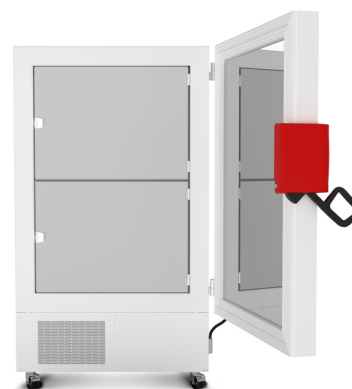

Model 700

IMPORTANT FEATURES

- Temperature range: -90 °C to -40 °C
- Lowest energy consumption in its class
- Energy Star certified (depending on version)
- Powerful cascade compressor cooling unit
- Environmentally friendly refrigerants R-290 and R-170
- Thermal insulation with unique long-life VIPs (vacuum insulation panels)
- Innovative door gasket concept reduces buildup of ice
- Ergonomic door handle
- Rust-proof interior made completely of stainless steel
- Detachable inner doors made of stainless steel
- 3 stainless-steel shelves providing flexible placement
- Water cooling (depending on version)
- Zero-voltage alarm contact for in-house alarm system
- Internal data logger, measured values can be read out in open format via USB
- Ethernet interface
- 3 access ports, Ø 28 mm, rear

ORDERING INFORMATION

Interior volume		Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
[cu.ft.]	[L]					
Model UF V 700						
26	736	230 V 1~ 50 Hz -10,0 A	CEE 7/7	Standard	UFV700-230V	9020-0348
				Water cooling	UFV700-230V-W	9020-0350
		115 V 1~ 60 Hz -15,0 A	NEMA 5-15	Standard	UFV700UL-120V	9020-0352
		208...230 V 2~ 60 Hz -10,0 A	NEMA 6-15	Standard	UFV700UL-240V	9020-0354


TECHNICAL DATA

Data				
Article Number	9020-0348	9020-0350	9020-0354	9020-0352
Designation	UFV700-230V	UFV700-230V-W	UFV700UL-240V	UFV700UL-120V
Option model	Standard	Water cooling	Standard	Standard
Performance Data Temperature				
Temperature range	-90...-40 °C	-90...-40 °C	-90...-40 °C	-90...-40 °C
Temperature uniformity at -80°C	2.5 ± K	2.5 ± K	2.5 ± K	2.5 ± K
Temperature fluctuation at -80°C	1.5 ± K	1.5 ± K	1.5 ± K	1.5 ± K
Cooling down time from 22°C to -80°C	450 min	450 min	450 min	450 min
Electrical data				
Rated Voltage	230 V	230 V	208...230 V	115 V
Power frequency	50 Hz	50 Hz	60 Hz	60 Hz
Nominal power	1,6 kW	1,6 kW	1,8 kW	1,4 kW
Unit fuse	10,0 A	10,0 A	10,0 A	15,0 A
Phase (Nominal voltage)	1~	1~	2~	1~
Dimensions and weights				
Interior volume	26 cu.ft.	26 cu.ft.	26 cu.ft.	26 cu.ft.
Net weight of the unit (empty)	635 lbs	635 lbs	635 lbs	635 lbs
Permitted load of individual compartment	110 lbs	110 lbs	110 lbs	110 lbs
Permitted load	441 lbs	441 lbs	441 lbs	441 lbs
Wall clearance back	3,94 in	3,94 in	3,94 in	3,94 in
Wall clearance sidewise	9,84 in	9,84 in	9,84 in	9,84 in
Housing dimensions not incl. fittings and connections				
Width net	43.7 in	43.7 in	43.7 in	43.7 in
Height net	77.4 in	77.4 in	77.4 in	77.4 in
Depth net	36.93 in	36.93 in	36.93 in	36.93 in
Internal Dimensions				
Interior width	35.04 in	35.04 in	35.04 in	35.04 in
Interior height	51.18 in	51.18 in	51.18 in	51.18 in

Interior depth	25.04 in	25.04 in	25.04 in	25.04 in
Inner doors	2	2	2	2
Unit doors	1	1	1	1
Environment-specific data				
Sound-pressure level	47 dB(A)	47 dB(A)	47 dB(A)	47 dB(A)
Average heat compensation at set value -80 °C and T _a = 21 °C	340 W	340 W	340 W	340 W
Energieverbrauch bei Sollwert -80 °C und T _u = 20 °C	8.1 kWh/d	8.1 kWh/d	8.1 kWh/d	8.1 kWh/d
Fixtures				
Number of shelves (std./max.)	3/11	3/11	3/11	3/11
Load-variants				
Freezer racks per compartment	6	6	6	6
Cryoboxes, 50 mm	528	528	528	528


All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS Incl. fittings and connections [mm]

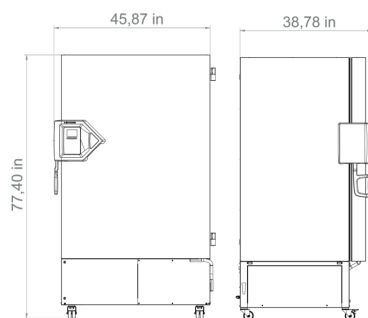

DIAGRAMS


Pull-up time during power failure

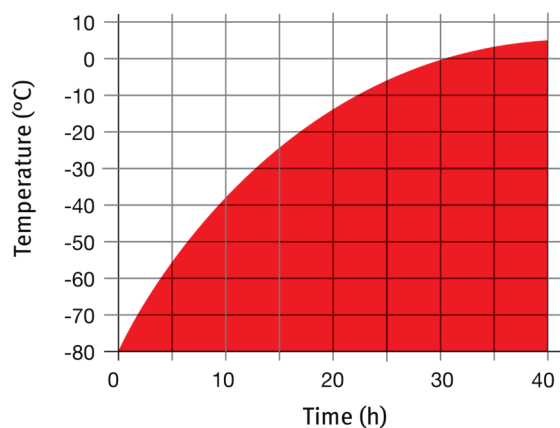

Ambient temperature 22 °C, no load

Heating up rate

Pull-down time

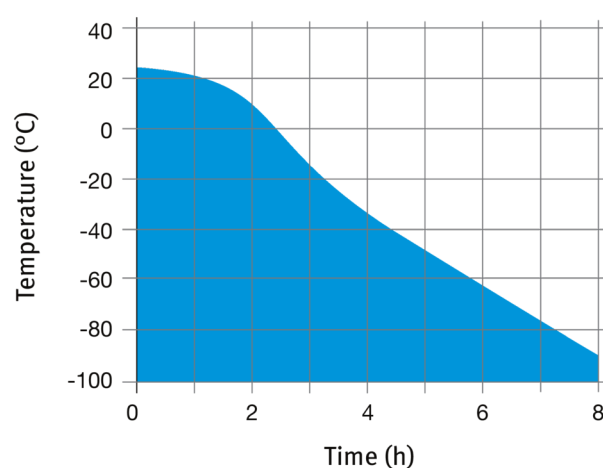

Ambient temperature 25 °C, no load

Cooling down rate

OPTIONS

Designation	Description	UF V 700	*	Art.-No.
Alarm system, battery-backed	for maintaining temperature display, alarm messages, and communication	•	-	8012-1910
Analog output 4-20 mA	for temperature values (output not adjustable)	•	-	8012-1898
Buck-boost converter	automatic voltage compensation	•	-	8012-1897
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	-	8012-1126
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	•	-	8012-0031
	temperature measurement incl. certificate, 20 measuring points at specified temperature	•	-	8012-0925
Compartment doors	Foam-filled compartment doors with gasket	•	-	8012-1894
CO ₂ backup cooling	factory installed (allow for increased depth of UF V: 1085 mm), includes the option of a battery-backed alarm system	•	-	8012-1936
Door access system	consisting of electromechanical door locking mechanism and electronic access control via NumPad	•	-	8012-1899
Door lock	lockable door handle	•	-	8012-1911
Shipping pallet	with unloading ramp for installing unit without forklift	•	-	8012-1907
Switch cover	Lockable main switch cover	•	-	8012-1980

ACCESSORIES

Designation	Description	UF V 700	*	Art.-No.
5 m power cable	C 13 refrigerator connection cable with BS 1363 British plug	•	-	5023-0285
	C 13 refrigerator connection cable with CEE 7/7 European plug	•	-	5023-0283
	C 13 refrigerator connection cable with type 13 Swiss plug	•	-	5023-0284
APT-COM™ 4 BASIC-Edition	for simple logging and documentation requirements with up to 5 networked units.			
	version 4, BASIC edition	•	-	9053-0039
APT-COM™ 4 GLP-Edition	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.			
	version 4, GLP edition	•	-	9053-0042
APT-COM™ 4 PROFESSIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.			
	version 4, PROFESSIONAL edition	•	-	9053-0040
Chamber rack, aluminum	for 5x4 standard boxes (50 mm), dimensions: H 278 x W 139 x D 564 mm			
	Set of 3	•	-	8012-2345
	Set of 4	•	-	8012-2346
	for 5x4 standard boxes (50 mm) incl. boxes, dimensions: H 278 x W 139 x D 564 mm			
	Set of 3	•	-	8012-2347
	Set of 4	•	-	8012-2348
	for 6x4 standard boxes (50 mm), dimensions: H 334 x W 139 x D 564 mm			
	Set of 3	•	-	8012-2357
	Set of 4	•	-	8012-2358
	for 6x4 standard boxes (50 mm) incl. boxes, dimensions: H 334 x W 139 x D 564 mm			
	Set of 3	•	-	8012-2359
	Set of 4	•	-	8012-2360
Chamber rack, stainless steel	for 5x4 standard boxes (50 mm), dimensions: H 278 x W 139 x D 561 mm			
	Set of 3	•	-	8012-2349
	Set of 4	•	-	8012-2350
	for 5x4 standard boxes (50 mm) incl. boxes, dimensions: H 278 x W 139 x D 561 mm			
	Set of 3	•	-	8012-2351

Designation	Description	UF V 700	*	Art.-No.
Chamber rack with drawer, stainless steel	Set of 4	•	-	8012-2352
	for 6x4 standard boxes (50 mm), dimensions: H 334 x W 139 x D 561 mm			
	Set of 3	•	-	8012-2361
	Set of 4	•	-	8012-2362
	for 6x4 standard boxes (50 mm) incl. boxes, dimensions: H 334 x W 139 x D 561 mm			
	Set of 3	•	-	8012-2363
	Set of 4	•	-	8012-2364
	for 5x4 standard boxes (50 mm), dimensions: H 278 x W 140 x D 562 mm			
	Set of 3	•	-	8012-2353
	Set of 4	•	-	8012-2354
	for 5x4 standard boxes (50 mm) incl. boxes, dimensions: H 278 x W 140 x D 562 mm			
	Set of 3	•	-	8012-2355
	Set of 4	•	-	8012-2356
	for 6x4 standard boxes (50 mm), dimensions: H 334 x W 140 x D 562 mm			
	Set of 3	•	-	8012-2365
	Set of 4	•	-	8012-2366
	for 6x4 standard boxes (50 mm) incl. boxes, dimensions: H 334 x W 140 x D 562 mm			
	Set of 3	•	-	8012-2367
	Set of 4	•	-	8012-2368
Cryoboxes	set consisting of 36 cardboard boxes with sample dividers	•	-	8012-2369
Deicing kit	consisting of de-icing tray, wiper and scraper	•	-	8012-1909
Grid shelf	For temperature range -40 °C to -60 °C, 3 provided in set, including shelf supports			
	stainless steel	•	-	8012-1904
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	-	8012-2250
Qualification documents IQ/OQ	IQ/OQ documents digital in PDF format – supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature. Excluded: Document execution is separate			
	temperature. Excluded: Document execution is separate	•	-	7057-0010

Designation	Description	UF V 700	*	Art.-No.
Qualification documents IQ/OQ/ PQ	IQ/OQ/PQ documents digital in PDF format – supporting documents for validation customized to customer requirements, PQ section added to the qualification folder IQ/OQ; supporting documents for validation consisting of: IQ/OQ checklists, comprehensive unit documentation covering: temperature values. Excluded: Document execution is separate			
	temperature values. Excluded: Document execution is separate	•	-	7057-0016
Shelf	Incl. shelf supports			
	stainless steel	•	-	8012-1902

SERVICES

Designation	Description	*	Art.-No.
Calibration services			
Temperature calibration	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate Travel to be quoted separate.	05	DL30-0190
Temperature measurement, 20 measuring points	Temperature measurement with 20 measuring points with a set value specified by the user, including certificate	05	DL30-0620
Installation services			
Unit installation	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) Travel to be quoted as separate.	05	DL10-0170
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming. Travel to be quoted as separate Limited to 3 hours.	05	DL10-0540
Maintenance services			
BRONZE 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate. 10% Discount on Spare Parts and Labor for any additional Services needed Wear + tear parts to be purchased or billed to customer. Travel to be quoted as separate.	05	DL20-0720
GOLD 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate of Temperature. 25% Discount on Spare Parts and Labor for any additional Services needed. One additional courtesy visit per year can be used for service visit or operator. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	05	DL20-0955
Preventive maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Checking a test temperature specified by the user in center of usable space without certificate. Wear + tear parts need to be purchased or billed to customer. Travel to be quoted as separate.	05	DL20-0603
SILVER 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate 15% Discount on Spare Parts and Labor for any additional Services needed One additional courtesy visit per year can be used for service visit or operator training. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	05	DL20-0855
Validation services			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder. Travel to be quoted separate.	05	DL41-0210

Designation	Description	*	Art.-No.
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder. Travel to be quoted separate.	05	DL44-0500

NOTES

- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany
TEL +49 7462 2005 0
info@binder-world.com
www.binder-world.com

BINDER Inc.

Bohemia, NY, USA
TEL +1 631 224 4340
usa@binder-world.com
www.binder-world.us

**BINDER Environmental Testing
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China
TEL +86 21 685 808 25
china@binder-world.com
www.binder-world.com