

Model WIC 3 | Walk-in-chambers for stability testing according to ICH (Q1A)


The walk-in-chambers offer a wide range of features and benefits to perfectly meet requirements in stability testing. Perfect test conditions are maintained in a wide range of temperature from 10°C to 50°C and relative humidity of 20% to 90% RH.

BENEFITS

- · Stable and uniform environment
- Fully complies with Pharma regulations (ICH)Equipment
- Longevity corrosion resistant
- · Ease of operation
- · Environmental friendly refrigerant with GWP value 547
- · Maintenance possible with the least possible disruption to the climate conditioning
- Perfectly meets requirements in Pharma industry
- · 3 years warranty



Model 3

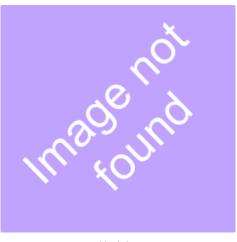

Model 3

IMPORTANT FEATURES

- Relative humidity accuracy of ± 2.5%
- Temperature accuracy of ± 1.5°C
- Duct work for proper distribution of air
- Highly insulated room panels with U value of only 0,21 $\mbox{W/m}^2\mbox{K}$
- LED Lighting with motion detector
- · High quality room elements, made in Germany
- Interior completely stainless steel

- · Exterior service accessibility
- Room Elements: Galvanized sheet steel, polyester coated, white, similar to RAL 9010, interior surfaces in stainless steel,
- 4 inch / 100 mm PU foam U value: $0,21W/m^2K$
- Floor: Floor overlay is chrome-nickel steel, pattern-rolled 5WL
- Surface load: 50.000 N/m²

ORDERING INFORMATION

Interior volume		Power supply - unit fuse	Plug*	Version	Model version	ArtNo.
[L]	[L]					
Model WIC 3						
24,516	24,516	230 V 1~ 50 Hz -20,0 A	fixed supply	Standard	WIC3-230V	9800-0013
		208230 V 1~ 60 Hz -20,0 A	fine decomply	Standard	WIC3-208V	9800-0007
			fixed supply	with conditioner on ceiling	WIC3C-208V 9800-0008	9800-0008




TECHNICAL DATA

TECHNICAL DATA					
Data					
Designation	WIC3-230V	WIC3-208V	WIC3C-208V		
Article Number	9800-0013	9800-0007	9800-0008		
Option model	Standard	Standard	with conditioner on ceiling		
Performance Data Climate					
Temperature range with humidity	1050 °C	1050 °C	1050 °C		
Humidity range	2090 % r.F.	2090 % r.F.	2090 % r.F.		
Temperature uniformity at 25°C and 60% RH	1.5 ± K	1.5 ± K	1.5 ± K		
Temperature uniformity at 40°C and 75% RH	1.5 ± K	1.5 ± K	1.5 ± K		
Temperature fluctuation at 25°C and 60% RH	0.5 ± K	0.5 ± K	0.5 ± K		
Temperature fluctuation at 40°C and 75% RH	0.5 ± K	0.5 ± K	0.5 ± K		
Electrical data					
Rated Voltage	230 V	208230 V	208230 V		
Power frequency	50 Hz	60 Hz	60 Hz		
Nominal power	2,7 kW	2,7 kW	2,7 kW		
Unit fuse	20,0 A	20,0 A	20,0 A		
Phase (Nominal voltage)	1~	1~	1~		
Measures					
Interior volume	24,516 L	24,516 L	24,516 L		
Net weight of the unit (empty)	1,186 kg	1,186 kg	1,186 kg		
Housing dimensions not incl. fittings and con	nections				
Width net	2,600 mm	2,600 mm	2,600 mm		
Height net	2,470 mm	2,470 mm	2,470 mm		
Depth net	4,700 mm	4,700 mm	4,700 mm		
Internal Dimensions					
Interior width	2,400 mm	2,400 mm	2,400 mm		
Interior height	2,270 mm	2,270 mm	2,270 mm		
Interior depth	4,500 mm	4,500 mm	4,500 mm		
Unit doors	1	1	1		




Environment-specific data			
Energy consumption at 40°C and 75% RH	1,150 Wh/h	1,150 Wh/h	1,150 Wh/h

DIMENSIONS Incl. fittings and connections [mm]

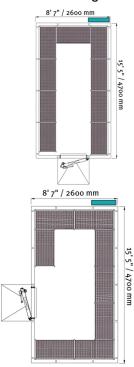

DIAGRAMS

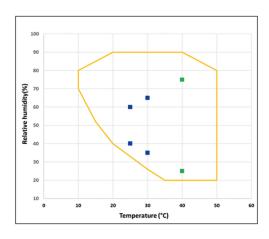

Temperature-environmental chart

OPTIONS

No data found.



ACCESSORIES

Designation	Description	WIC 3	*	ArtNo.
Qualification documents	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature and humidity values			
	Hard copy inside folder	•	-	7007-0002
Description all MIG mandals	Suitable for all WIC models	•	-	8012-2155
Ramp for all WIC models		•	-	8012-2239
GL IS ALL	suitable for WIC 3, WIC 3 C	•	-	8012-2158
Shelf set large		•	-	8012-2242

SERVICES

Designation	Description	*	ArtNo.			
	1 access port with silicone plug, 30 mm - position of your choice	01	DL10-1000			
	1 access port with silicone plug, 50 mm - your choice of position	01	DL10-1050			
	2 access ports with silicone plugs, 30 mm - your choice of position	01	DL10-1010			
	2 access ports with silicone plugs, 50 mm - your choice of position	01	DL10-1060			
Access to the effect of the ef	3 access ports with silicone plugs, 30 mm - your choice of position	01	DL10-1020			
Access port with silicone plug	3 access ports with silicone plugs, 50 mm - your choice of position	01	DL10-1070			
	4 access ports with silicone plugs, 30 mm - your choice of position	01	DL10-1030			
	4 access ports with silicone plugs, 50 mm - your choice of position	01	DL10-1080			
	5 access ports with silicone plugs, 30 mm - your choice of position	01	DL10-1040			
	5 access ports with silicone plugs, 50 mm - your choice of position	01	DL10-1090			
Calibration services						
Temperature measurement, 13 temperature measuring points and 1 humidity measuring point	Temperature measurement with 13 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate	-	DL30-0113			
Installation services	Installation services					
Installation and Comissioning	Installation and Comissioning	-	DL10-0800			
Maintenance services						
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	-	DL20-1000			
Validation services						



Designation	Description	*	ArtNo.
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder. Travel to be quoted separate.	-	DL40-0200



NOTES

O1 Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.

BINDER GmbH

Tuttlingen, Germany
TEL +49 7462 2005 0
FAX +49 7462 2005 100
info@binder-world.com
www.binder-world.com

BINDER Inc.

Bohemia, NY, USA
TEL +1 631 224 4340
FAX +1 631 224 4354
usa@binder-world.com
www.binder-world.us

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd.

Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com