

## Model BF 400 | Standard-Incubators with forced convection

The BINDER incubator of the BF series with mechanical convection is suitable for all gentle incubation applications, particularly under a full load and large batch throughputs. This incubator is extremely homogeneous and has quick recovery times.

### BENEFITS

- Safe and reproducible incubation even under high batch throughputs in long-term operation
- Disinfection routine at 100 °C
- Stable atmosphere by inner glass door



Model 400

### MAIN FEATURES

- Temperature range: room temperature plus 5 °C to 100 °C
- APT.line™ preheating chamber technology
- Adjustable fan speed
- Adjustable exhaust air flap
- Controller with timer functions
- Inner door made of tempered safety glass
- 2 chrome-plated racks
- Class 3.1 independent temperature safety device (DIN 12880) with visual temperature alarm
- Computer interface: RS 422

### ORDERING INFORMATION


| Interior volume [L] | Voltage              | Option model | Version      | Art.-No.  |
|---------------------|----------------------|--------------|--------------|-----------|
| <b>Model BF 400</b> |                      |              |              |           |
| 400                 | 230 V 1~ ph 50/60 Hz | Standard     | BF400-230V   | 9010-0241 |
|                     | 120 V 1~ ph 60 Hz    | Standard     | BF400UL-120V | 9010-0242 |

## TECHNICAL DATA

| Description  | BF400-230V <sup>1</sup> | BF400UL-120V <sup>1</sup> |
|--|-------------------------|---------------------------|
| Article Number   | 9010-0241               | 9010-0242                 |
| <b>Performance Data Temperature</b>                      |                         |                           |
| Temperature range 5 °C above ambient temperature to [°C] | 100                     | 100                       |
| Temperature variation at 37 °C [± K]                     | 0.4                     | 0.4                       |
| Temperature fluctuation at 37 °C [± K]                   | 0.1                     | 0.1                       |
| Heating-up time to 37 °C [min]                           | 12                      | 12                        |
| Recovery time after 30 seconds door open at 37 °C [min]  | 5                       | 5                         |
| <b>Electrical data</b>                                   |                         |                           |
| Rated Voltage [V]  | 230                     | 120                       |
| Power frequency [Hz]                                     | 50/60                   | 60                        |
| Nominal power [kW]                                       | 0.85                    | 0.85                      |
| Unit fuse [A]  | 12.5                    | 12.5                      |
| Phase (Nominal voltage)                                  | 1~                      | 1~                        |
| <b>Outer dimensions</b>                                  |                         |                           |
| Width net [mm]   | 1235                    | 1235                      |
| Height net [mm]  | 1025                    | 1025                      |
| Depth net [mm]   | 765                     | 765                       |
| Wall clearance back [mm]                                 | 160                     | 160                       |
| Wall clearance sidewise [mm]                             | 100                     | 100                       |
| <b>Doors</b>   |                         |                           |
| Unit doors   | 2                       | 2                         |
| <b>Internal Dimensions</b>                               |                         |                           |
| Width [mm]   | 1000                    | 1000                      |
| Height [mm]  | 800                     | 800                       |
| Depth [mm]   | 500                     | 500                       |
| <b>Measures</b>  |                         |                           |
| Interior volume [L]                                      | 400                     | 400                       |
| Net weight of the unit (empty) [kg]                      | 145                     | 145                       |
| permitted load [kg]                                      | 90                      | 90                        |
| Load per rack [kg]                                       | 35                      | 35                        |
| <b>Environment-specific data</b>                         |                         |                           |
| Energy consumption at 37 °C [Wh/h]                       | 56                      | 56                        |
| <b>Fixtures</b>  |                         |                           |
| Number of shelves (std./max.)                            | 2/9                     | 2/9                       |

1 All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

## DIMENSIONS incl. fittings and connections [mm]

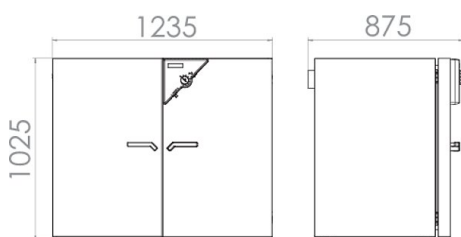

## OPTIONS AND ACCESSORIES

| Designation                          | Description  | *         | Art.-No.  |
|--------------------------------------|--|-----------|-----------|
| Access port with silicone plug       | back   |           |           |
|                                      | 10 mm  | 01        | 8012-1281 |
|                                      | 30 mm  | 01        | 8012-1291 |
|                                      | 50 mm  | 01        | 8012-1293 |
|                                      | left   |           |           |
|                                      | 10 mm  | 01        | 8012-1274 |
|                                      | 30 mm  | 01        | 8012-1156 |
|                                      | 50 mm  | 01        | 8012-1164 |
|                                      | right  |           |           |
|                                      | 10 mm  | 01        | 8012-1267 |
|                                      | 30 mm  | 01        | 8012-1153 |
|                                      | 50 mm  | 01        | 8012-1162 |
|                                      | top  |           |           |
|                                      | 10 mm  | 01        | 8012-1261 |
|                                      | 30 mm  | 01        | 8012-1150 |
|                                      | 50 mm  | 01        | 8012-1159 |
| 100 mm                               | 01, 10   | 8012-1166 |           |
| Alarm function for overheating       | Switchable acoustic alarm, with adjustable limit value on the independent temperature safety device  |           |           |
|                                      | 120 V option model   | -         | 8012-1302 |
|                                      | 230 V option model   | -         | 8012-1298 |
| Analog output 4-20 mA                | for temperature values (output not adjustable)   | -         | 8012-0482 |
| APT-COM™ 4 BASIC-Edition             | for simple logging and documentation requirements with up to 5 networked units.  |           |           |
|                                      | version 4, BASIC edition   | 19        | 9053-0039 |
| APT-COM™ 4 GLP-Edition               | for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.                                    |           |           |
|                                      | version 4, GLP edition   | 19        | 9053-0042 |
| APT-COM™ 4 PROFESSIONAL-Edition      | convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.   |           |           |
|                                      | version 4, PROFESSIONAL edition  | 19        | 9053-0040 |
| Calibration certificate, expanded    | for temperature; for extending the measurement in center of chamber to include another test temperature  | -         | 8012-1111 |
|                                      | for temperature, measurement in center of chamber at specified temperature   | -         | 8012-1130 |
| Calibration certificate, temperature | temperature measurement incl. certificate and 27 measuring points at specified temperature   | -         | 8012-1590 |
|                                      | temperature measurement incl. certificate, 15- 18 measuring points at specified temperature  | -         | 8012-1569 |
|                                      | temperature measurement incl. certificate, 9 measuring points at specified temperature   | -         | 8012-1548 |
|                                      |  | -         |           |
| Data Logger Kit                      | T 220: For continuous temperature logging from -90 °C to 220 °C. The kit includes 1 data logger, Pt 100 sensor with 2 m extension cable and 1 magnetic fixture for mounting to the BINDER unit | 19        | 8012-0715 |
| Data Logger Software                 | LOG ANALYZE software kit, configuration and evaluation software for all BINDER Data Logger Kits (incl. USB data cable)   | 19        | 8012-0821 |
| Door lock                            | lockable door handle   | -         | 8012-1252 |
| Inner chamber, reinforced            | including 2 reinforced racks   | -         | 8012-1732 |
| Interface converter                  | RS 422 cable set and RS 232 / RS 422 interface converter for connection to 10-way plug distributor or for connection to a unit   |           |           |
|                                      | 120 V, 60 Hz option model  | -         | 8012-0557 |
|                                      | 230 V, 50/60 Hz option model   | -         | 8012-0556 |
|                                      | RS 422 cable set and RS 422 / Ethernet interface converter for connection to 10-way plug distributor or for connection to a unit   |           |           |
|                                      | 120 V, 60 Hz option model  | -         | 8012-0405 |
|                                      | 230 V, 50/60 Hz option model   | -         | 8012-0380 |
| Interface converter, USB / RS422     | RS 422 cable set and RS 422 / USB interface converter for connection to 10-way plug distributor or for connection to a unit (USB-powered converter)  | -         | 8012-0665 |

\* Notes &gt; See last page

| Designation                    | Description  | *  | Art.-No.  |
|--------------------------------|--|----|-----------|
| Interior socket                | Waterproof, switched flush-mounting box (CA3GD) for nominal voltage in unit interior, with cover and corresponding plug (nominal voltage; max. 500 W; max. 90 °C; protection class IP67) | 07 | 8012-1170 |
| pH-neutral detergent           | concentrated, for gentle remove of residual contaminants; 1 kg   | -  | 1002-0016 |
| Pt 100 temperature sensor      | additional flexible Pt 100, interior, with external connection incl. LEMO connector (3-pin)  | -  | 8012-1172 |
| Rack                           | chrome plated  | -  | 6004-0005 |
|                                | stainless steel  | -  | 6004-0011 |
| RS 422 interface               | modular plug distributor for 10 RS 422 interfaces  | -  | 8012-0295 |
| RS 422 interface, cable (15 m) | RS 422 connection cable (15 m) between plug distributor and RS 422 interface   | -  | 5023-0036 |
| RS 422 interface, cable (50 m) | RS 422 extension cable (50 m) between interface converter and unit or RS 232 / RS 422 plug distributor   | -  | 5023-0117 |
| Rubber pads                    | set anti-slip feet   | -  | 8012-1887 |
| Shelf, perforated              | Stainless steel  | -  | 6004-0032 |
| Table on castors               | stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm   | -  | 9051-0019 |

\* Notes › See last page

## SERVICES

| Designation                                 | Description   | *              | Art.-No.  |
|---|---|----------------|-----------|
| <b>Installation services</b>                |   |                |           |
| Instruction                                 | unit function instructions for operation and programming of the controller  | 18             | DL10-0500 |
| <b>Preventive maintenance</b>               |   |                |           |
| Preventive maintenance                      | Executive of equipment inspection according to maintenance plan   | 14, 18         | DL20-0200 |
| <b>Calibration services</b>                 |   |                |           |
| Calibration temperature incl. certificate   | including certificate, one measuring point in center of chamber at specified temperature  | 14, 16, 17, 18 | DL30-0101 |
| Temperature measurement 18 measuring points | including certificate, 18 measuring points at specified temperature   | 14, 16, 17, 18 | DL30-0118 |
| Temperature measurement 27 measuring points | including certificate, 27 measuring points at specified temperature   | 14, 16, 17, 18 | DL30-0127 |
| Temperature measurement 9 measuring points  | including certificate, 9 measuring points at specified temperature  | 14, 16, 17, 18 | DL30-0109 |
| <b>Validation services</b>                  |   |                |           |
| Execution of IQ/OQ                          | in accordance with qualification folder   | 15, 18, 20     | DL40-0100 |
| Execution of IQ/OQ/PQ                       | in accordance with customer's requirement, price: on request  | 15, 18         | DL44-0500 |
| Qualification folder IQ/OQ                  | supporting documents for validation performed by customer, consisting of: IQ/OQ checklists, unit schematics, QM certificate in accordance with ISO 9001     | 15, 18, 20     | 8012-0871 |
| Qualification folder IQ/OQ/PQ               | supporting documents for validation performed by customer, in accordance with customer's requirement, extension of Qualification folder IQ/OQ by chapter PQ | 15, 18         | 8012-0959 |
| <b>Warranty service</b>                     |   |                |           |
| Extension of the warranty from 2 to 3 years | beginning with the date of delivery, wearing parts are not included   | -              | DL01-7041 |
| Extension of the warranty from 2 to 5 years | beginning with the date of delivery, wearing parts are not included   | -              | DL01-7042 |

\* Notes › See last page

**NOTES**

- 01 Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 02 UL mark is not granted when this option is used.
- 03 Heat resistant only to max. 200 °C.
- 04 Only available on units rated for 230 V.
- 06 Heating-up time may increase as a result of the lower heat conductivity.
- 07 The additional heat input may affect the temperature behavior.
- 09 Not in conjunction with the optional access port, door with window and interior lighting.
- 10 Not available on 23-liter units.
- 11 Not available on 23- or 53-liter units.
- 12 Only available on units rated for 230 V or 400 V.
- 13 Installation and connections take place at unit location; transport within the company only upon consultation.
- 14 We recommend a BINDER service contract to cover unit inspections, calibrations and validations.
- 15 OQ according to Yellow Paper = completed factory validation documentation of all OQ checklists.
- 16 Sensor calibration is performed in an accredited calibration laboratory.
- 17 Calibration is performed according to the BINDER factory standard.
- 18 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).
- 19 For additional accessories, refer to the Process documentation chapter.
- 20 When ordering IQ/OQ qualification folders and associated IQ/OQ execution on one order, we offer a 15 % discount for the item of the IQ/OQ folder.  
When ordering IQ/OQ/PQ qualification folders and associated IQ/OQ/PQ execution on one order, we offer a 15 % discount for the item of the IQ/OQ/PQ folder.

BINDER GmbH  
Tuttlingen, Germany  
TEL +49 7462 2005 0  
FAX +49 7462 2005 100  
info@binder-world.com  
www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.  
Kowloon, Hong Kong, P.R. China  
TEL +852 39070500  
FAX +852 39070507  
asia@binder-world.com  
www.binder-world.com

Representative Office for CIS  
Moscow, Russia  
TEL +7 495 988 15 16  
FAX +7 495 988 15 17  
russia@binder-world.com  
www.binder-world.ru

BINDER Environmental Testing  
Equipment (Shanghai) Co., Ltd.  
Shangha, P.R. China  
TEL +86 21 685 808 25  
FAX +86 21 685 808 29  
china@binder-world.com  
www.binder-world.com

BINDER Inc.  
Bohemia, NY, USA  
TEL +1 631 224 4340  
FAX +1 631 224 4354  
usa@binder-world.com  
www.binder-world.us