


Measuring device



Individual extras for your simulation chamber

ADVANTAGES FOR CUSTOMERS

- APT-COM™ recording of non-BINDER units. The sensors in this measuring device are incorporated into the chamber. There is an Ethernet communication interface.
- Additional monitoring of the chamber interior. Direct monitoring of the test specimen is also available, depending on the sensor technology required.
- 4 sensors (temperature, humidity, CO₂, O, etc.) can be displayed and recorded at the same time.
- Set values can be entered, which means that if the adjustable tolerance band is exceeded or not reached, there will be a visual and acoustic alarm.
- The device is attached and removed purely magnetically, so no structural measures are required on the chamber.

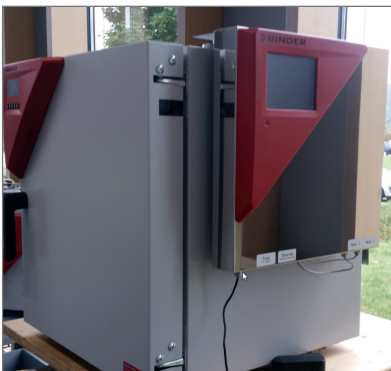

BINDER INDIVIDUAL inquiries
go2binder.com/en-BINDER-INDIVIDUAL